116949

- CARBON, TOTAL ORGANIC AND INORGANIC


Carbon may exist in sediment and water samples as either inorganic or organic compounds. Inorganic carbon is present as carbonates, bicarbonates, and possibly free carbon dioxide. Specific types of compounds that are considered to be included in the organic carbon fraction are nonvolatile organic compounds (sugars), volatile organic compounds (mercaptans), partially volatile compounds (oils), and particulate carbonaceous materials (cellulose).^{1,2[#]}

The basis of the method is the catalytic or chemical oxidation of carbon in carbon-containing compounds to carbon dioxide followed by the quantification of the carbon dioxide produced. Alternately, the carbon may be reduced to methane and appropriately quantified. It follows, then, that the distinction between inorganic carbon and organic carbon is the method of sample pretreatment. There are presently two procedures for defining this separation. One method is based on sample treatment with a strong acid. Analysis of an untreated sample is a measure of total carbon while analysis of the acid-treated fraction is a measure of organic carbon. Inorganic carbon is calculated by subtraction. The second method of separation is based on differential thermal combustion with organic compounds being converted to carbon dioxide at $500^{\circ}C$ to $650^{\circ}C^{3+4}$ and inorganic carbon being converted to carbon dioxide at $950^{\circ}C$ to $1300^{\circ}C$.

Sample Handling and Storage

Flowcharts for the handling of samples intended for organic carbon and inorganic carbon analysis are presented in Figure 3-6 and Figure 3-7. Water and sediment samples to be analyzed for inorganic carbon may be stored in glass or plastic containers. There is no effective preservative because of the carbon dioxide reserve in the atmosphere. The only precaution that can be taken for inorganic

* References for this procedure can be found on page 3-76. AR301362

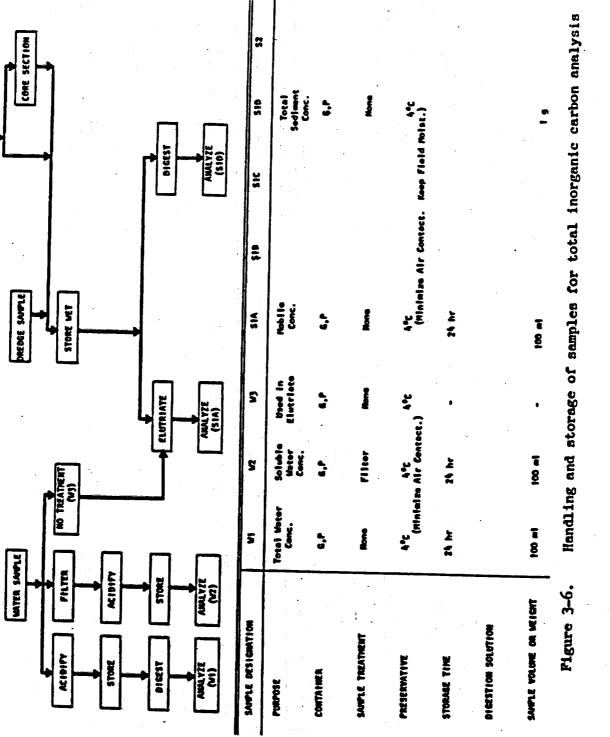

3

ł

CORE SAMPLE


• : :

AR301363

.

. .

C

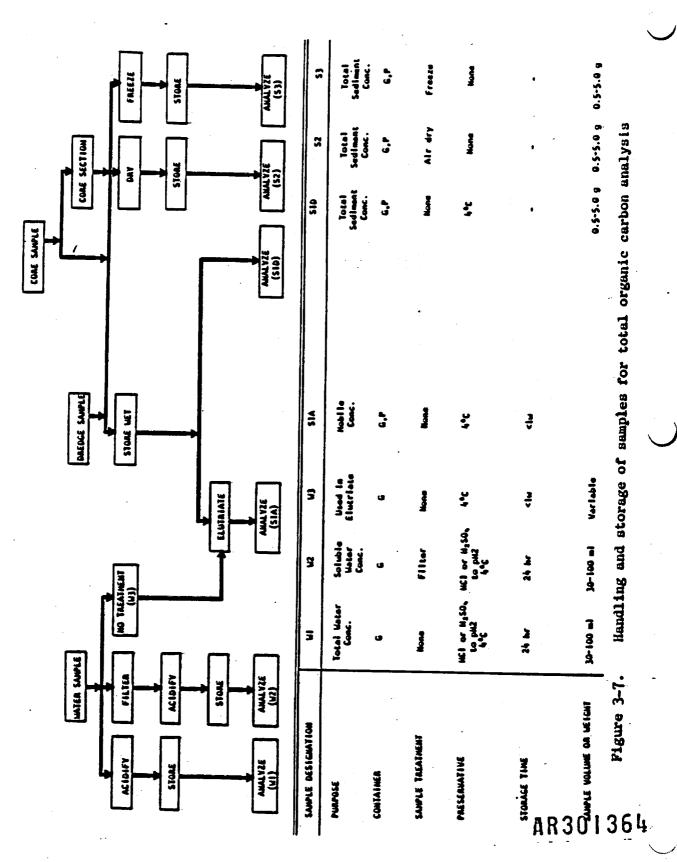

carbon is to completely fill the sample container at the time of sampling (exclude all air bubbles), tightly seal the container, and complete the analysis immediately (Figure 3-6).

Water samples for organic carbon analysis should be stored in glass containers unless substitute containers have been shown not to affect total organic carbon (TOC) analyses. Samples should be processed as soon as possible (within 24 hr if possible) to minimize change due to chemical or biological oxidation. Atmospheric uptake of carbon dioxide is less critical since it would be evolved when the sample is acidified prior to analysis. Sediment samples for organic carbon analysis may be stored in either plastic or glass containers (Figure 3-7). Air drying of sediments (S2) may lead to low TOC results due to oxidation or volatilization. Therefore, moist storage (SID) or frozen storage (S3) would be the preferred method of storage. If samples are frozen, excessive temperatures should not be used to thaw the samples.

AR301365

: 3-68

Procedure for Water Samples (W1, W2, S1A)

Method 1: Infrared Analysis^{5,7}

Apparatus

Sample homogenizer such as a Waring blender or ultrasonic blender

Magnetic stirrer

Hypodermic syringe

Total carbon analyzer, either a single channel or a dual channel instrument (Dov-Beckman Carbonaceous Analyzer Model No. 915, Dohrmann Envirotech DC-50 carbon analyzer, Oceanography International Total Carbon Analyzer, Leco, or equivalent)

Reagents

Distilled water: the distilled water used in the preparation of standards and dilution of samples should be of the highest quality in order to have a small blank.

Organic carbon, stock solution, 1000 mg/f C: dissolve 2.125 g anhydrous potassium biphthalate, KHC: H.O., in distilled water and dilute to 1 ¹/₂ in a volumetric flask.

Organic carbon, standard solutions: prepare standard solution by dilution of the stock solution as required.

Inorganic carbon, stock solution, 1000 mg/2: dissolve 3.500 g sodium bicarbonate, NaHCO3, and 4.418 g sodium carbonate, Na2CO3, in distilled water in a 1-2 volumetric flask and make up to the mark.

Inorganic carbon, standard solution: prepare standards from the stock solution as required.

Packing for total carbon tube: dissolve 20 g cobalt nitrate, $CO(NO_3)_2$ · $6H_2O_1$, in 50 ml distilled water. Add this solution to 15 g long-fiber asbestos in a porcelain evaporating dish. Mix and evaporate to dryness on a steam bath. Place the dish in a muffle furnace and bring to 950°C. After 1 to 2 hr at this temperature, remove the dish and allow to cool. Break up any large lumps and mix adequately but not excessively. With the combustion tube held in a vertical position, taper joint up, put about 1/2 in. of untreated asbestos in the tube first, then transfer in small amounts, approximately 1 g of catalyst into the tube with forceps or tweezers. As it is added, tap or push the material gently with a 1/4-in. glass rod. Do not force the packing. The weight of the rod itself is sufficient to compress the material. When completed, the length of the packing should be about 30 of 355 6 cm. Test the packed tube by measuring the flow rate of

gas through it at room temperature, and then at 750°C. The rate should not drop more than 20 percent.

 $\mathbf{h} \in \mathbb{R}^{+}$

Packing for carbonate tube (dual channel instrument): place a small wad of quartz wool or asbestos near the exit end of the carbonate evolution tube. From the entrance end add 6 to 12 mesh quartz chips, allowing these to collect against the wad to a length of 10 cm. Pour an excess of 85 percent phosphoric acid, H₃PO₄, into the tube while holding it vertically and allow the excess to drain out.

Nitrogen gas, carbon dioxide free.

. ja

Procedure

Turn on the infrared analyzer, recorder, and tube furnaces, setting the total carbon furnace at 950°C and the carbonate furnace at 175°C. Allow sufficient warm-up time for stable, drift-free operation; about 2 hr is required. If used daily, the analyzer can be left on continuously. Adjust the oxygen flow rate to 80 to 100 ml/min through the total carbon tube. With other instruments, follow manufacturer's directions to warm up the instrument.

Immediately prior to carrying out calibrations or analyses, inject several portions of the appropriate standard into the tube to be used, until constant readings are obtained. The actual injection technique is as follows: rinse the syringe several times with the solution to be analyzed, fill, and adjust the volume to be pipeted. Wipe off the excess with soft paper tissue, taking care that no lint adheres to the needle. Remove the plug from the syringe holder, insert the sample syringe, and inject the sample into the combustion tube with a single, rapid movement of the thumb. Leave the syringe in the holder until the flow rate returns to normal, then replace it with the plug.

Successively introduce a convenient sized aliquot (20 to $50 \ \mu$) of each organic carbon standard and a blank into the total carbon tube and record peak heights. Between injections allow the recorder pen to return to its baseline. When a dual channel instrument is used, the standardization procedure must be repeated using carbonate standards to calibrate the low temperature channel.

Thoroughly mix the sample. Inject a convenient sized aliquot (20 to 50 μ) of the sample into the total carbon tube and AR301367

record the peak height. This result is a measure of the organic carbon concentration and the inorganic carbon concentration of the sample.

Thoroughly mix the sample using a Waring blender or an ultrasonic homogenizer. Transfer 10 to 15 ml of sample to a 30-ml beaker and acidify with concentrated HCl to a pH of 2 or less. Purge the sample with carbon dioxide free nitrogen gas for 5 to 10 min. Plastic tubing should not be used during the purging process unless it has been previously shown that it will not add organic carbon to the sample.

Mix the acidified sample on a magnetic stirrer. While stirring, withdraw a subsample from the beaker using a hypodermic needle with a 150-µm opening. Inject the sample into the carbon analyzer to be used and record the peak height. This result is a measure of the organic carbon concentration of the sample.

Using either clear or filtered water samples, analytical precision will approach 1 to 2 percent or 1 to 2 mg/l carbon, whichever is greater. Analytical precision for unfiltered water samples will increase to 5 to 10 percent because of the difficulty associated with sampling particulate matter and the fact that the needle opening of the syringe limits the maximum size of the particles that can be included in the sample.

Calculations

<u>Dual-channel instrument.</u> Prepare calibration curves derived from the peak heights obtained with the standard total carbon and inorganic carbon solutions.

Determine the concentration of total carbon and inorganic carbon in the sample by comparing sample peak heights with the calibration curves.

Determine the concentration of total inorganic carbon in the sample by subtracting the organic carbon value from the total carbon value.

Single-channel instrument. Prepare a calibration curve derived from the peak heights obtained with the standard total carbon solutions. Determine the total carbon concentration in the sample by comparing the peak height of the first sample injection with 363363

calibration curve. Determine the organic carbon concentration in the sample by comparing the peak height of the second sample injection with the calibration curve. Inorganic carbon concentrations are calculated by subtracting the organic carbon concentration from the total carbon concentration.

<u>*</u>.

Procedures for Sediment Samples (S1D, S3)

Method 1: Sample Ignition

Apparatus

Induction furnace such as the Leco WR-12, Dohrmann DC-50, Coleman CH analyzer, or Perkin Elmer 240 elemental analyzer

Combustion boats

Microbalance

Desiccator

Reagents

10 percent hydrochloric acid: mix 100 ml concentrated HCl with 900 ml distilled water.

Copper oxide fines.

Benzoic acid.

Procedure

Dry at 70°C and grind the sediment sample.

Weigh a combustion boat and record the weight. Place 0.2 to 0.5 g homogenized sediment in the combustion boat and reweigh. Combustion boats should not be handled with the bare hand during this process.

If total carbon or inorganic carbon is to be determined, Cupric oxide fines may be added to the sample to assist in combustion. Combust the sample in an induction furnace. Record the result as total carbon.

If organic carbon is to be determined, treat a known weight of dried sediment with several drops of 10 percent HCL. Wait until the effervescing is completed and add more acid. Continue this process until the incremental addition of acid causes no further effervescence. Do not add too much acid at one time as this may cause loss of sample due to frothing.

Dry the sample at 70°C and place in a desiccator. Add Cupric oxide fines, combust the sample in an induction furnace, and record the result as organic carbon.

AR301370

Calculations

The carbon content of the sample can be calculated as:

SC = weight of tube (after-before) x 27.29 sample weight

Derivation of factor:

27.29 = <u>12.011 (molecular weight carbon)</u> × 100%

When the total sample results are used, the result is percent carbon in the sample. When acid-treated samples are used, the result is percent organic carbon. Inorganic carbon is calculated as total carbon minus organic carbon.

Method 2: Differential Combustion4,5

Apparatus

Sargent programmed microcombustion apparatus or equivalent Microbalance

Procedure

.....

-1642 - 1642 - Air dry the sediment sample. Using a mortar and pestle, grind the sample to pass a 100-mesh screen.

Combust a known weight of sediment at a programmed heating rate of 300° to 950° C in 10 min and then maintain 950° C for 20 min. Trap the CO₂ in ascarite and record the weight as total carbon. A sample size should be selected that will produce 25 to 50 mg CO₂.

Weigh a second portion of the dried sediment. Combust this sample at a programmed rate of 300° to 650° C in 10 min and maintain 650° C for 20 min. Trap the CO₂ in ascarite and record the weight as organic carbon.

Calculations

The total carbon concentration, C_t , of the sample (in mg/g) is calculated as follows:

$$C_{t} = \frac{(x_{t}) \left(\frac{12}{44}\right)}{(g)}$$

where

xt = weight of CO₂ evolved at 950°C, mg

g = weight of sample combusted, g

The organic carbon, C_0 , concentration of the sample (in $\pm g/g$) is calculated as follows:

$$C_{0} = \frac{(x_{0}) (\frac{12}{64})}{(g)}$$

where

xo = weight of CO₂ evolved at 650°C, mg
g = weight of sample combusted, g

Inorganic carbon, C_{I} , (in mg/g) is calculated as: $C_{I} = C_{t} - C_{o}$

Method 3: Wet Combustion***

A third method has been used for carbon in sediments. This is based on the oxidation of the sample with dichromate and back titration of the sample with ferrous ammonium sulfate. References are provided for the procedure but details are not given. The procedure is similar to the chemical oxygen demand test which is not specific for carbon. The vet combustion method is a redox procedure and any reduced chemicals in the sediment samples (ferrous iron, manganous manganese, sulfide) will react with the dichromate. Therefore, this procedure is not recommended unless other instrumentation is not available.

References

- 1. U. S. Environmental Protection Agency. "Manual of Methods for Chemical Analysis of Water and Wastes." Methods Development and Quality Assurance Research Laboratory, National Environmental Research Center; Cincinnati, Ohio. 298 p. (1974).
- 2. U. S. Environmental Protection Agency. "Methods for Chemical Analysis of Water and Wastes." Environmental Monitoring and Support Laboratory, Office of Research and Development, IPA; Cincinnati, Ohio (1979).
- Giovannini, G., Poggio, G., and Sequi, P. "Use of an Automatic CHN Analyzer to Determine Organic and Inorganic Carbon in Soils." Unpublished Report, Laboratory of Soil Chemistry, via Corridoni, Pisa, Italy. 9 p. (1975).
- 4. Konrad, J. G., Chesters, G., and Keeney, D. R. "Determination of Organic- and Carbonate-Carbon in Freshwater Lake Sediments by a Microcombustion Procedure." J. Thermal Analysis 2:199-208 (1970).
- Kemp, A. L. W. "Organic Matter in the Sediments of Lakes Ontario and Erie." Proc. 12th Conference Great Lakes Research 12:237-249 (1969).
- 6. Environment Canada. "Analytical Methods Manual," Inland Waters Directorate, Water Quality Branch; Ottawa, Canada (1974).
- American Public Health Association. <u>Standard Methods for the Examination of Water and Wastevater</u>. APHA; New York, New York. 1193 p. (1976).
- Gaudette, H. E., Flight, W. R., Toner, L., and Folger, D. W.
 "An Inexpensive Titration Method for the Determination of Organic Carbon in Recent Sediments." J. Sed. Fetrology 44:249-253 (1974).

METHOD 9045

SOIL pH

1.0 SCOPE AND APPLICATION

1.1 Method 9045 is an electrometric procedure which has been approved for measuring pH in calcareous and noncalcareous soils.

2.0 SUMMARY OF METHOD

2.1 The soil sample is mixed either with Type II water or with a calcium chloride solution (see Section 5.0), depending on whether the soil is considered calcareous or moncalcareous. The pH of the solution is then measured with a pH meter.

3.0 INTERFERENCES

3.1 Samples with very low or very high pH may give incorrect readings on the meter. For samples with a true pH of >10, the measured pH may be incorrectly low. This error can be minimized by using a low-sodium-error electrode. Strong acid solutions, with a true pH of $\langle 1, may give incorrectly$ high pH measurements.

3.2 Temperature fluctuations will cause measurement errors.

3.3 Errors will occur when the electrodes become coated. If an electrode becomes coated with an oily material that will not rinse free, the electrode can either (1) be cleaned with an ultrasonic bath, or (2) be washed with detergent, rinsed several times with water, placed in 1:10 HCl so that the lower third of the electrode is submerged, and then thoroughly rinsed with water.

4.0 APPARATUS AND MATERIALS

4.1 pH Meter with means for temperature compensation.

4.2 Electrodes:

4.2.1 Calonel electrode.

4.2.2 Glass electrode.

4.2.3 A combination electrode can be employed fastead of calomel or glass.

4.5 Beaker: 50-mL.

9045 - 1

AR301374 Revision <u>0</u> Date September 1986 4.6 Volumetric flask: 2-Liter.

4.7 <u>Volumetric flask</u>: 1-Liter.

5.0 REAGENTS

5.1 <u>ASTM Type II water</u> (ASTM D1193): Water should be monitored for impurities.

5.2 <u>Primary standard buffer salts</u> are available from the National Bureau of Standards (NBS) and should be used in situations where extreme accuracy is necessary. Preparation of reference solutions from these salts requires some special precautions and handling, such as low-conductivity dilution water, drying ovens, and carbon-dioxide-free purge gas. These solutions should be replaced at least once each month.

5.3 <u>Secondary standard buffers</u> may be prepared from NBS salts or purchased as solutions from commercial vendors. These commercially available solutions, which have been validated by comparison with NBS standards, are recommended for routine use.

5.4 <u>Stock calcium chloride solution</u> (CaCl₂), 3.6 M: Dissolve 1059 g of CaCl₂·2H₂O in Type II water in a 2-liter volumetric flask. Cool the solution, dilute it to volume with Type II water, and mix it well. Dilute 20 mL of this solution to 1 liter with Type II water in a volumetric flask and standardize it by titrating a 25-mL aliquot of the diluted solution with standard 0.1 N AgNO₃, using 1 mL of 5% K₂CrO₄ as the indicator.

5.5 <u>Calcium chloride</u> (CaCl₂), 0.01 H: Dilute 50 mL of stock 3.6 M CaCl₂ to 18 liters with Type II water. If the pH of this solution is not between 5 and 6.5, adjust the pH by adding a little Ca(OH)₂ or HCl. As a check on the preparation of this solution, measure its electrical conductivity. The specific conductivity should be 2.32 ± 0.08 mmho per cm at 25°C.

6.D SAMPLE COLLECTION, PRESERVATION, AND HANDLING

6.1 All samples must be collected using a sampling plan that addresses the considerations discussed in Chapter Nine of this manual.

5.2 Samples should be analyzed as soon as possible.

7.D PROCEDURE

7.1 Calibration:

7.1.1 Because of the wide variety of pH meters and accessories, detailed operating procedures cannot be incorporated into this method. Each analyst must be acquainted with the operation of each system and familiar with all instrument functions. Special attention to care of the electrodes is recommended.

9045 - 2

AR301375 Revision 0 Date September 1986 7.1.2 Each instrument/electrode system must be calibrated at a minimum of two points that bracket the expected pH of the samples and are approximately three pH units or more apart. Repeat adjustments on successive portions of the two buffer solutions until readings are within 0.05 pH units of the buffer solution value.

7.2 Sample preparation and pH measurement of noncalcareous soils:

7.2.1 To 20 g of soil in a 50-mL beaker, add 20 mL of Type II water and stir the suspension several times during the next 30 min.

7.2.2 Let the soil suspension stand for about 1 hr to allow most of the suspended clay to settle out from the suspension.

7.2.3 Adjust the electrodes in the clamps of the electrode holder so that, upon lowering the electrodes into the beaker, the glass electrode will be immersed just deep enough into the clear supernatant solution to establish a good electrical contact through the ground-glass joint or the fiber-capillary hole. Insert the electrodes into the sample solution in this manner. For combination electrodes, immerse just below the suspension.

7.2.4 If the sample temperature differs by more than 2°C from the buffer solution, the measured pH values must be corrected.

7.2.5 Report the results as "soil pH measured in water."

7.3 Sample preparation and pH measurement of calcareous soils:

7.3.1 To 10 g of soil in a 50-mL beaker, add 20 mL of 0.01 M CaCl₂ (Step 5.5) solution and stir the suspension several times during the next 30 min.

7.3.2 Let the soil suspension stand for about 30 min to allow most of the suspended clay to settle out from the suspension.

7.3.3 Adjust the electrodes in the clamps of the electrode holder so that, upon lowering the electrodes into the beaker, the glass electrode will be immersed well into the partly settled suspension and the calomel electrode will be immersed just deep enough into the clear supernatant solution to establish a good electrical contact through the ground-glass joint or the fiber-capillary hole. Insert the electrode into the sample solution in this manner.

7.3.4 If the sample temperature differs by more than 2°C from the buffer solution, the measured pH values must be corrected.

7.3.5 Report the results as "soil pH measured in 0.01 H CaCl₂".

9045 - 3

Revision 0 Date September 1986

8.0 QUALITY CONTROL

8.1 Duplicate samples and check standards should be analyzed routinely.

8.2 Electrodes must be thoroughly rinsed between samples.

9.0 METHOD PERFORMANCE


9.1 No data provided.

10.0 REFERENCES

11.11

10.1 None required.


AR301377 Revision 0 Date September 1986

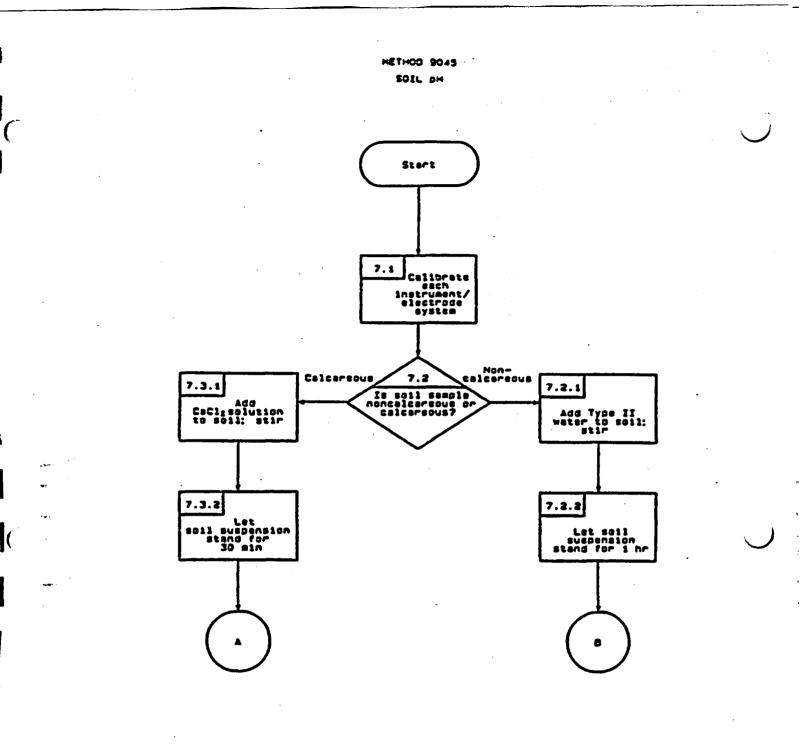


0

Date September 1986

Revision

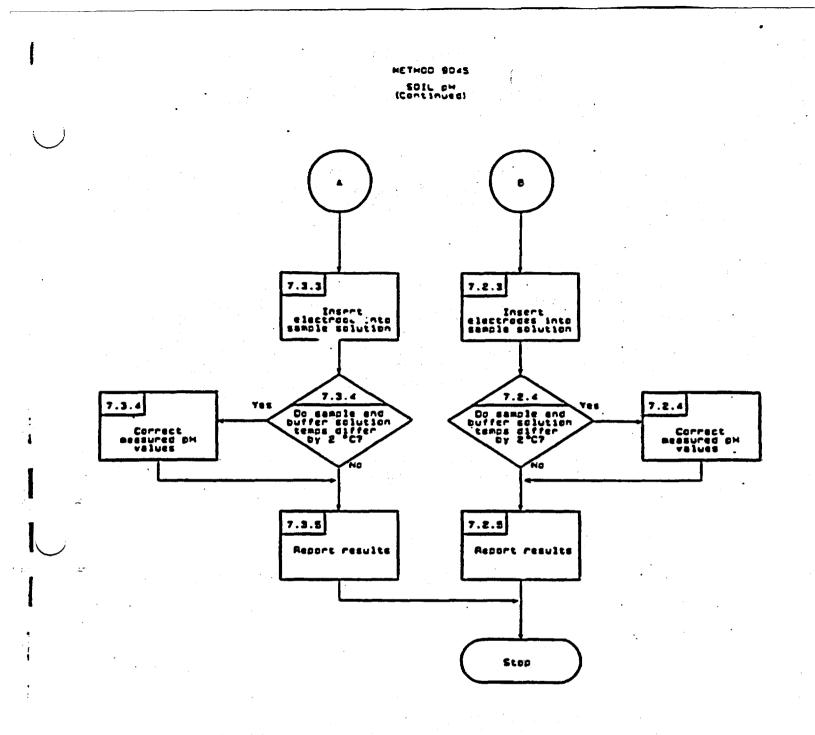

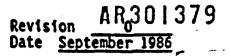

METHOD 9081

CATION-EXCHANGE CAPACITY OF SOILS (SODIUM ACETATE)

1.0 SCOPE AND APPLICATION

1.1 Method 9081 is applicable to most soils, including calcareous and noncalcareous soils. The method of cation-exchange capacity by summation (Chapman, 1965, p. 900; see Paragraph 10.1) should be employed for distinctly acid soils.

2.0 SUMMARY OF METHOD

2.1 The soil sample is mixed with an excess of sodium acetate solution, resulting in an exchange of the added sodium cations for the matrix cations. Subsequently, the sample is washed with isopropyl alcohol. An ammonium acetate solution is then added, which replaces the adsorbed sodium with ammonium. The concentration of displaced sodium is then determined by atomic absorption, emission spectroscopy, or an equivalent means.

3.0 INTERFERENCES

é.

3.1 Interferences can occur during analysis of the extract for sodium content. Thoroughly investigate the chosen analytical method for potential interferences.

4.0 APPARATUS AND MATERIALS

4.1 Centrifuge_tube and stopper: 50-mL, round-bottom, narrow neck.

4.2 <u>Mechanical shaker</u>.

4.3 Volumetric flask: 100-mL.

5.0 REAGENTS

5.1 Sodium acetate (NaOAc), 1.0 N: Dissolve 136 g of NaC₂H₂O₂·3H₂O in water and dilute it to 1,000 mL. The pH of this solution should be 8.2. If needed, add a few drops of acetic acid or NaOH solution to bring the reaction of the solution to pH 8.2.

5.2 Annonium acetate (NH40Ac), 1 N: Dilute 114 mL of glacial acetic acid (99.5%) with water to a volume of approximately 1 liter. Then add 138 mL of concentrated ammonium hydroxide (NH40H) and add water to obtain a volume of about 1,980 mL. Check the pH of the resulting solution, add more NH40H, as needed, to obtain a pH of 7, and dilute the solution to a volume of 2 liters with water.

9081 - 1

AR301380

Revision <u>O</u> Date September 1986

5.3 Isopropyl alcohol: 99%.

• 6.0 SAMPLE COLLECTION, PRESERVATION, AND HANDLING

6.1 All samples must be collected using a sampling plan that addresses the considerations discussed in Chapter Nine of this manual.

7.0 PROCEDURE

> ، ملی ... بر این ... بر این ...

7.1 Weigh 4 g of medium- or fine-textured soil or 6 g of coarse-textured soil and transfer the sample to a 50-mL, round-bottom, narrow-neck centrifuge tube. (A fine soil has >50% of the particles <0.074 mm, medium soil has >50% >0.425 mm, while a coarse soil has more than 50\% of its particles >2 mm.

7.2 Add 33 mL of 1.0 N NaOAc solution, stopper the tube, shake it in a mechanical shaker for 5 min, and centrifuge it until the supernatant liquid is clear.

7.3 Decant the liquid, and repeat Paragraph 7.2 three more times.

7.4 Add 33 mL of 99% isopropyl alcohol, stopper the tube, shake it in a mechanical shaker for 5 min, and centrifuge it until the supernatant liquid is clear.

7.5 Repeat the procedure described in Paragraph 7.4 two more times.

7.6 Add 33 mL of NH4OAc solution, stopper the tube, shake it in a mechanical shaker for 5 min, and centrifuge it until the supernatant liquid is clear. Decant the washing into a 100-mL volumetric flask.

7.7 Repeat the procedure described in Paragraph 7.6 two more times.

7.8 Dilute the combined washing to the 100-mL mark with ammonium acetate solution and determine the sodium concentration by atomic absorption, emission spectroscopy, or an equivalent method.

8.0 QUALITY CONTROL

8.1 All quality control data should be maintained and available for easy reference or inspection.

8.2 Employ a minimum of one blank per sample batch to determine if contamination or any memory effects are occurring.

8.3 Materials of known cation-exchange capacity must be routinely analyzed.

9081 - 2

Revision A<u>R301381</u> Date September 1986

9.0 METHOD PERFORMANCE


9.1 No data provided.

10.0 REFERENCES

í

10.1 This method is based on Chapman, H.D., "Cation-exchange Capacity," pp. 891-900, in C.A. Black (ed.), Method of Soil Analysis, Part 2: Chemical and Microbiological Properties, Am. Soc. Agron., Madison, Wisconsin (1965).



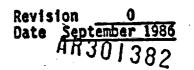


ŗ

NETHOD 908 1-24 CATION-EXCHANGE CAPACITY OF BOILS (BODIUN ACETATE)

-

K)

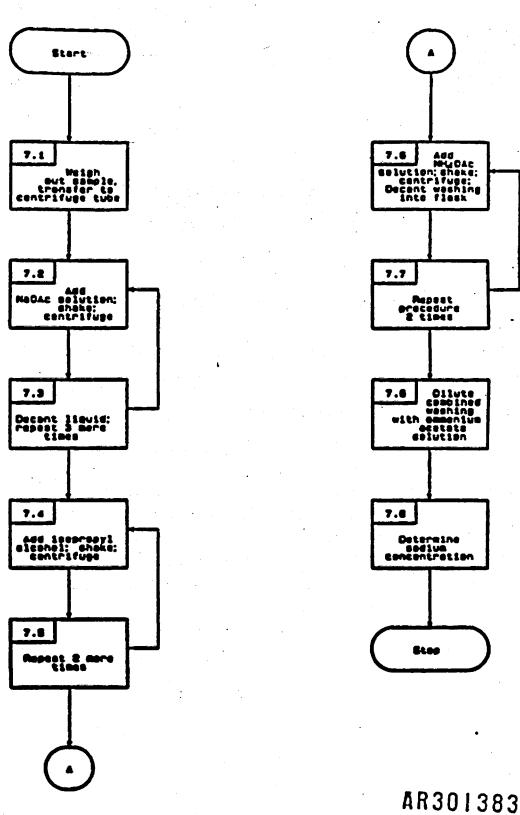

Revision 0 Date September 1986

....

---- - -

REM III SAS	No.	
553	09	03

AR301384

Ebasco Services Incorporated REM III Program Zone Program Management Office 2000 15th Street North Arlington, VA 22201 Phone: (703) 558-7555

REM III TEAM LABORATORY SPECIAL ANALYTICAL SERVICES Regional Office Request

	Regional Transmittal
A.	Regional Office: <u>NUS-Pittsburch</u> and Site Name: Croydon TCE
8.	Regional Laboratory Service Coordinator: Greg Zimmerman
Ç.	Telephone Mumber: (412) 788-1080
-152 D.	Date of Request: September 22, 1988

Please provide below a description of your request for Special Analytical Services under the REM III Program. In order to most efficiently obtain capability for your request, please address the following considerations, if applicable. Incomplete or erroneous information may result in delay in processing your request. Provide response on additional sheets or attach supplementary information as needed. 1. General description of analytical service requested: Analysis of 4 Soil samples for BIU, Percent Ash, and TOC. 2. Definition and number of work units involved (specify whether whole samples or fractions; whether organics or inorganics; whether aqueous or soil and sediments; or whether low, medium, or high concentration: 4 low concentration soil samples. 3. Purpose of analysis (specify whether Remedial or Enforcement Action, RCRA, NPDES, etc.): RI/FS 4. Data Quality Objectives (specify whether screening analysis, engineering support or confirmational analysis): Engineering Support DOO Level III. 5. Estimated date(s) of collection: November 7 - 11, 1988

- 6. Estimated date(s) and method of shipment: November 7 11, 1988. Samples will be shipped by overnight air carrier. These dates are tentative and are dependent on the project remaining on schedule.
- 7. Approximate number of days results are required after receipt of samples by laboratory: 40 days
- 8. Analytical protocol required (attach copy if protocol other than currently approved under this contract or CLP program): BTU - ASTM D 3286-85 All methods are attached. Percent Ash - ASTM D 2974-84 TOC - page 3-73 from "Procedures for Handling and Chemical Analysis of Sediment and Water Samples", EPA, May 1981.
- 9. Special technical instructions (if outside protocol requirements, specify compound names, detection limits, etc.):_____
- 10. Analytical results required (if known, specify format for data sheets, QA/QC reports, chain-of-custody, documentation, etc.). If not completed, format of results will be left to discretion of the laboratory: Data sheets, chain-of-custody forms, all calculations, and SAS Request Form.
- 11. Other instructions (use additional sheets or attach supplementary information, as needed):__

-AR301386

12. Name of sampling/shipping contact: Greg Zimmerman

AR301387

Phone: (412) 788-1080

I. DATA REQUIREMENTS

	mg/l	PRECISION DESIRED
PARAMETER	DETECTION LIMIT	(+ % or Conc.)
BIU	Not Specified	Not Specified
Percent Ash	Not Specified	Not Specified
TOC	Not Specified	Not Specified
	- <u></u>	<u></u>
	·····	
- <u></u>		


II. QUALITY CONTROL REQUIREMENTS

SAMPLES	CONTROL	PRECISION DESIRED
CONTROL REQUIRED	FREQUENCY OF SAMPLES	LIMITS (+ or Conc.)
None Specified		·
	- <u></u>	
		-

III. ACTION REQUIRED IF LIMITS ARE EXCEEDED:

None Specified

TAR301388

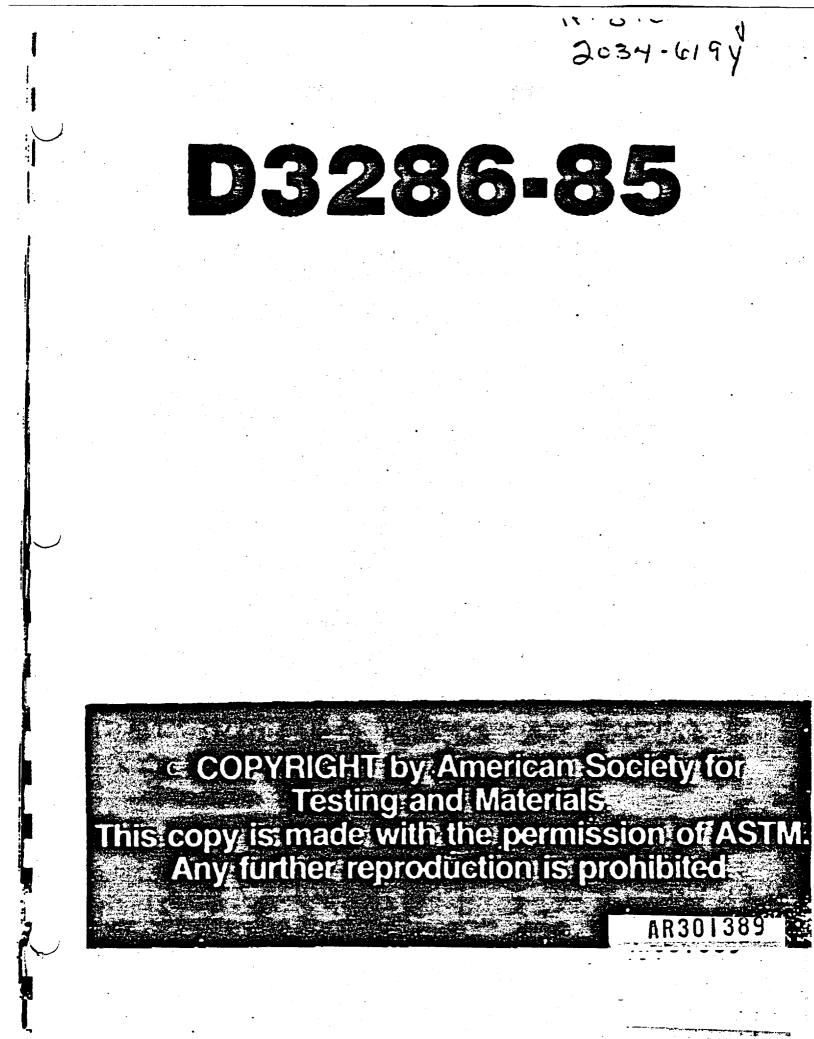

4 001007

M Designation: D 3286 - 65

Standard Test Method for GROSS CALORIFIC VALUE OF COAL AND COKE BY THE ISOTHERMAL-JACKET BOMB CALORIMETER'

This standard is issued under the fixed designation D 3214, the number immediately following the designation indicates the year of singular adoption or in the case of revision, the year of last revision. A number in parenthetes indicates the year of last reupproval A superscript epidon to indicates an editional change since the last revision or reapproval.

1. Scope-

1.1 This test method covers the determination of the gross calorific value of coal and coke by the isothermal-jacket bomb calorimeter.

1.2 The values stated in SI units are to be regarded as the standard.

1.3 This standard may involve hatardous indterials, operations, and equipment. This standard does not purport to address all of the safety problemis associated with its use. It is the responsibility of whoever uses this standard to consult and establish appropriate safety and health practices and determine the applicability of regulatory limitations prior to use Specific precautionary statements are given in Section 9.

2. Applicable Documents-

2

1

÷.

- 2.1 ASTM Standards:
- D121 Definitions of Terms Relating to Coal and Coke²
- D 346 Method of Collection and Preparation of Coke Samples for Laboratory Analysis¹
- D1193 Specification for Reagent Water³
- D 2013 Method of Preparing Coal Samples for Analysis³
- D.1173 Test Method for Moisture in the Analysis Sample of Coal and Coke³
- D3177 Test Methods for Total Sulfur in the Analysis Sample of Coul and Coke³
- D.1180 Method for Calculating Coal and Coke Analyses from As-Determined to Different Baser?
- D4239 Test Method for Sulfur in the Analysis Sample of Coal and Coke Using High Temperature Tube Furnave Combustion Methods³
- E. I. Specification for ASTM Thermometers⁴
- E 144 Recommended Practice for Safe Use of Oxygen Combustion Bombs'

3. Terminology

3.1 Definitions:

3.1.1 calurific value—the heat produced by combustion of a unit quantity of a substance under specified conditions. It is expressed in this test method in British thermal units per pound (Btu/Ib). Calorific value may also be expressed in calories per gram (cal/g) or in the International System of Units (SI), joules per gram (J/g), when required. The unit equivalents are given in Table 1.

3.1.2 gross calorific value (gross heat of comhustion at constant volume), Q. (gross)—see Definitions D 121.

3.1.3 net calorific value (net heut of combustion at constant pressure). Q. (net)—see Definitions D 121.

3.1.4 calorimeter—as used in this test method, consists of the bomb and its contents, the calonmeter vessel (bucket) with stierer, the water in which the bomb is immersed, and the portions of the thermometer and the ignition leads within the valorimeter vessel.

3.2 Descriptions of Terms Specific to this Standard

3.2.1 corrected temperature race—the temperature of the calorimeter caused by the process that occurs inside the bomb; that is the observed temperature change corrected for various effects as noted in 10.4.1.

⁴ This test method is under the particle time of ASTAT conanitee 184 on Cital and Code and is the direct responsibility of Subsymmittee 1019 (1) on Methods of Anglesis

Subsymmittice EEPS 21 on Methods of Anglesis Current edition approved Aug. 23, 1985. Published Deuber 1984. Originally published in D.3286 - 73. Last previous edition-15 USB - 84.

Inning Bod of 1871 Standards, Not 05-05

* Innust Head of INTAT Ston Aude. Vet 11 01

Innind Back of 181 18 Standards, Vil 14 11

Innual Book of INEXE Standards And 14102

NOTE 1—Temperature is measured in either degrees Celsus or degrees Fahrenheit. Thermometer corrections should be applied. Temperatures may be recorded in ohms or other arbitrary units instead of degrees. Consistent units must be used in standardization and the actual calorific value determination. If arbitrary units other than degrees Celsus or Fahrenheit are used, the temperature interval over which all tests are made must not vary so much that an error greater than 0 KH1°C would be caused.

3.2.2 energy equivalent, heat capacity, or water equivalent—the energy required to raise the temperature of the calorimeter one arbitrary unit. This is the quantity that, when multiplied by the corrected temperature rise, then adjusted for extraneous heat effects, and divided by the mass of the sample, gives the gross calorific value.

Nots 2—Energy units for quantities listed throughout this test method are such that the number of energy units per gram of sample corresponds exactly to the number of Bu's per pound of sample. For brevny these are referred to as Bu's. The actual energies are smaller than those stated by the ratio of the number of pounds per gram (1/433.39). The energy equivalent of the calonmeter has the units (Blu/Ib) times (g/deg). Conversion to other units is discussed in Appendix X2. Time is expressed in minutes. Mass is expressed in grams.

4. Summary of Method

t an c

4.1 Calorific value is determined in this test method by burning a weighed sample, in oxygen, in a calibrated isothermal-jucket bomb calorimeter under controlled conditions. The calorimeter is standardized by burning benzoic acid. The calonific value of the sample is computed from temperature observations made before, during, and after combustion, and making proper allowances for heat contributed by other processes, and for thermometer and thermochemical corrections.

Note 3—Oxidation after sampling of susceptible tow-rank cost or lignite may result in a reduction of catorific value. Unnecessary exposure of the sample to air from the sime of sampling or delay in analysis shall be avoided.

5. Significance and Use -

5.1 The gross calorific value is used to compute the total calorific content of the quantity of coal represented by the sample for payment purposes, provided the buyer and the seller mutually agree upon this.

5.2 The gross calorific value is used in computing the calorific value versus suffur content to determine if the coal meets the regulatory requirements for industrial fuels.

5.3 The gross calorific value may be used for

D 3286 -

evaluating the effectiveness of beneficiation processes, or for research purposes.

6. Apparatus and Facilities

6.1 Test Spare—A room or area free fromdrafts and that can be kept at a reasonably uniform temperature for the time required for the determination. The apparatus shall be shielded from direct sunlight and radiation from other heat sources. Thermostatic control of room temperature and controlled relative humidity are desirable.

6.2 Combustion Bomb, shall be constructed of materials that are not affected by the combussion process or products sufficiently to introduce measurable heat input of alteration of end products. The bomb must be designed so that all liquid combustion products can be completely recovered by washing the inner surfaces. There must be no gas leakage during a ten. The bomb must be capable of withstanding a hydrostatic pressure test of 20 MPs (3000 psig) at room temperature without stressing any part beyond fts elaric limit.

6.3 Balance-A is boratory balance having the capability to weigh the sample to the nearest 0.0001 g. The balance should be checked periodically to determine its accuracy.

6.4 Calorimeter Vessel, made of metal with a tarnish-resistant coating, and with all outer surfaces highly polished. Its size shall be such that the bomb will be completely immersed in water when the calorimeter is assembled. It shall have a device for stirring the water thoroughly and at a uniform rate but with minimum heat input. Continuous stirring for 10 min shall not raise the calorimeter temperature more than 0.01°C (0.02°F) starting with identical temperatures in the calorimeter, room, and jacket. The immersed portion of the stirrer shall be coupled to the outside through a material of low-heat conductivity.

6.5 Jacket, a double-walled, air, or water-filled jacket. The calorimeter shall be insulated from the jacket by an air space. The sides, top, and bottom of the calorimeter vessel shall be apprusimately 10 mm from the inner wall of the jacket to minimize convection currents. Mechanical supports for the calorimeter vessel shall provide as little thermal conduction as possible. The jacket shall be capable of maintaining the temperature constant to within ±0.1°C (±0.2°F) of room temperature at a calorimeter temperature 2°C (4°F) below, and 2°C (4°F) or more above

AR301391

newn temperature throughout the determination. If a water-filled packet is used, it shall have a device for stirring the water at a uniform rate with minimum heat input:

6.6. Thermometers, used to measure temperatures in the calorimeter and jacket shall be of any of the following types or combination thereof.

6.5.1 Liquid-in-Gluss Thermometers, conforming to the requirements for ASTM Thermometers 56C, 56F, 116C, or 117C, 25 prescribed in Specification E1. The thermometers shall be tested for accuracy against a known standard (preferably by the National Bureau of Standards). For thermometers 56C and 56F, the calibration should be at intervals no larger than 2.0°C or 2.5°F over the entire graduated scale. The maximum difference in correction between any two test points shall be no more than 0.02°C or 0.05°F. For thermometers 116C and 117C. the calibration should be at intervals no larger than 0.5°C over the entire calibrated range. The maximum difference in correction between any two test points shall not be more than 0.02°C.

6.6.2 Berkmann Enformatial Thermometers, (glass-enclosed scale, adjustable) having a range of approximately 6°C in 0.01°C subdivisions reading upward and conforming to the requirements for Thermometer 115C, as prescribed in Specification E.1. may be used, hach of these thermometers shall be tested for accuracy against a known standard (preferably by the National Bureau of Standards) at intervals no larger than 1°C over the entire graduated scale. The maximum difference in the correction between any two test points shall not be more than 0.02°C.

6.6.3 Other Thermometers, of an accuracy cyuul to or better than 0.001°C, such as platinum resistance or linear thermistor thermometers are satisfactory and may be used if properly calibrated. A Wheatstone bridge and galvanometer capable of measuring resistance to 0.0001 fl are necessary for use with 25 fl platinum resistance thermometers.

6.7 Thermometer Accessive — A magnifier is required for reading liquid-in-glass thermometers to one tenth of the smallest scale division. This shall have a lens and holder designed so as to introduce no significant errors due to parallas.

6.8 Sample Holder—An open crucible of platinum, quartz, or acceptable base-metal alloy. Base-metal alloy crucibles are acceptable if, after a few preliminary firings, the weight does not change significantly between tests.

6.9 Ignition Wire—The ignition wire shall be 100 mm of 0.16-mm diameter (No. 34 B & S gage) nickel-chromlum alloy. (Chromel C) alloy. or iron wire. Platinum wire or palladium. 0.10mm diameter (No. 38 B & S gage) may be used. provided constant ignition energy is supplied. The length or mass of Ignition wire shall remain constant for all calibrations and calorific value determinations.

6.10 Ignition Circuit. for ignition purposes shall provide a 6 to 16 V alternating or direct current to the ignition wire. An ammeter or pilot light is required in the circuit to indicate when current is flowing. A step-down transformer, connected to an alternating current lighting circuit, or batteries, may be used.

6.11 Burer, used for the acid titration, shall have 0.1-mL divisions.

6.12 Automated Controller and Temperature Measuring Accessories, may be used.

7. Resgests-

7.1 Reagent Water, conforming to Type II of Specification D 1193, shall be used for preparation of reagents and washing of the bomb interior.

7.2 Purity of Reagents-Reagent grade chemicals shall be used in all tests. Unless otherwise indicated, it is intended that all reagents shall conform to the specification of the Committee on Analytical Reagents of the American Chemical Society, where such specifications are available.⁶ Other grades may be used, provided it is first ascertained that the reagent is of sufficiently high purity to permit its use without lessening the accuracy of the determination.

7.3 Benzue Acid (C,H₅COOH), shall be the National Bureau of Standards benzoie acid. The crystals shall be pelleted before use. Commercially prepared pellets may be used provided they are made from National Bureau of Standards benzoic acid. The value of heat of combustion of benzoic acid for use in the calibration calculations shall be in accordance with the value listed in the National Bureau of Standards certificate issued with the standard.

74 Method Orange, Method Red, or Method

**Respirat Chemicals: American Chemical Survey Specifications," Am Chemical Sar, Washington, DC. I or supprisons on the training of respiration local by the American Chemical Survey, we "Respirat Chemicals and Standards," In Joseph Roam, D. Yan Mostella Ch., Inc., New York, NV, and the "United States Pharmacopera."

9 W |

Purple Indicator, may be used to titrate the acid formed during comhustion. The indicator used shall be the same for both calibration and calorific value determinations.

7.5 Oxygen manufactured from liquid air, guaranteed to be greater than 99.5 % pure, should be used. Oxygen made by the electrolytic process may contain a small amount of hydrogen rendering it unfit without purification.

7.6 Suchum Curbonate Standard Solution, sodium carbonate (Na₂CO₃) should be dried for 24 h at 105°C. Dissolve 20.9 g in water and dilute to 1 L. One millilitre of this solution is equivalent to 10.0 Btu in the nitric acid (HNO₃) titration.

8. Sample -

1.

-

5.

た

-

8.1 The sample shall be the material pulverized to pass No. 60 (250- μ m) sieve, prepared in accordance with either Method D 346 for coke. or Method D 2013 for coal.

8.2 A separate portion of the sample should be analyzed simultaneously for moisture content in accordance with Method D 2013 and Test Method D 3173, so that calculations to other bases can be made.

8.3 Sulfur analysis shall be made in accordance with Test Methods D 3177.

9. Safety Precautions

9.1 The following precautions are recommended for safe calorimeter operation. Additional precautions are given in Recommended. Practice E 144. Also consult the calorimeter equipment manufacturer's installation and operating instructions before using the calorimeter.

9.2 The mass of coal or coke sample and the pressure of the oxygen admitted to the bomb must not exceed the bomb manufacturer's recummendations.

9.3 Bomb parts should be inspected carefully after each use. Threads of the main closure should be checked frequently for wear. Cracked or significantly worn parts should be replaced. The bomh should be returned to the manufacturer occasionally for inspection and possible proof firing.

9.4 The oxygen supply cylinder should be equipped with an approved type of sofety device, such as a reducing valve, in addition to the needle valve and pressure gage used in regulating the oxygen feed to the bomb. Valves, gages, and gashets must meet industry safety code. Suitable reducing valves and adaptors for 3 to 4-MPa (300 to 500-psi) discharge pressure are obtainable from commercial sources of compressed gas equipment. The pressure gage shall be checked periodically for accuracy.

9.5 During ignition of a sample, the operator must not permit any portion of her/his body to extend over the calorimeter.

9.6 When combustion aids are employed, extreme caution must be exercised not to exceed the bomb manufacturer's recommendations and to avoid damage to the bomb. Do not fire loose fluffy material, such as unpelleted benzoic acid, unless thoroughly mixed with the coal sample.

9.7 Do not fire the bomb if the homb has been dropped or turned over after loading, or there is evidence of gas leakage when the bomb is submerged in the calorimeter water.

9.8 For manually operated calorimeters, the ignition circult switch shall be of the momentary double-contact type, normally open, except when held closed by the operator. The switch should be depressed only long enough to fire the charge.

10. Standardization

10.1 The calorimeter is standardized by combustion of benzole acid.

10.2 Determine the energy equivalent of the culorimeter for a specific temperature rise as the average of a series of ten individual runs made over a period of not less than J days nor more than S days. To be acceptable, the standard deviation of the series shall be 6.5 Btu/°C (3.6 Btu/ °F) or less (see Table 2). For this purpose, any individual test may be discarded only if there is evidence indicating incomplete combustion. If this limitation is not met, investigate for the source of the problem, currect it, then repeat the entire series to obtain a standard deviation within the acceptable limits.

10.3 Prixedure:

10.3.1 Regulate the weights of the peliets of benzole acid in each series to yield approximately the same temperature rise as that obtained with the coals tested in the same laboratory. The usual range of masses is from 0.9 to 1.3 g. Weigh the pellet to the nearest 0.0001 g in the sumple holder in which it is to be burned, and record the weight as the mass.

10.3.2 Rinse the bomh, invert to drain, and have undried. Add 1.0 mL of water to the homb

prior to assembly for a determination.

10.3.3 Connect a measured length of ignitionwire to the ignition ferminals, with enough slack to allow the ignition wire to maintain contact with the sample.

10.3.4 Assemble the bomb and charge it with oxygen to a consistent pressure between 2 to 3. MPa (20 to 30 atm). This pressure must remain the same for each calibration and each calorificvalue determination. Admit the oxygen slowly into the bomb so as not to blow powdered material from the sample holder. If the pressure exceeds the specified pressure, do not proceed with the combustion Instead, detach the filling connection and exhaust the bomb in the usual manner, and discard the sample.

10.3.5 Fill the calorimeter vessel (bucket) with the measured (or weighed) quantity of water adjusted from 1.0 to 2.0°C (2.0 to 4.0°F) below the jacket temperature, but not lower than 20°C (68°F) (Note 4). Use the same mass of water in each test weighed to ±0.5 g. For 2000-mL calorimeters, the proper quantity can be obtained by use of a volumetric flash calibrated to deliver 2000 ± 0.5 mL. As the density of water varies with temperature, make suitable corrections if the water temperature varies from the temperature at which the flash was calibrated. Place the assembled homb in the calorimeter vessel. Check that no oxygen hubbles are leaking from the homb. Place the calorimeter vessel in the jacket: connect the electrodest place the stirrers, there mometers, and cover in position. Start the stircerts) and continue to operate it (them) throughout the determination. Examine the thermometers for liquid separation and correct any separation before proceeding. The starting temperature should be within ±0.5°C (±0.9°F) of that used in the analysis of coal or coke samples.

Next 4— The initial temperature adjustment will ensure a final temperature slightly above that of the jacket for 2000-mil, calorimeters. Some operators prefer a lower initial temperature so that the final temperature is slightly, below that of the jacket. This procedure is also statifactory. Whichever procedure is used, the same procedure should be used in all tests, including standardization. A small heater may be built into the calorimeter so that the desired starting temperature can be easily attained.

Note 5-Check all liquid-in-glass thermometers at least daily for any defects, for example, cracked glass, etc.

10.3.6 Allow 5 min for attainment of equilibrium; then record the calorimeter temperature at t-min intervals for 5 min or until the rate of change has been nearly constant for 5 min. Use a magnifier when using ASTM Bomb Calonmeter Thermumeters 560 or 56E and estimate all

a magnifier when using ASTM Bomb Calonmeter Thermometers 56C or 56F, and estimate all of the readings (except those during the rapidrise period) to the nearest 0.002°C or 0.005°F. Estimate ASTM Thermometers 115C, 116C, or 117C readings to 0.001°C, and 25 12 resistance thermometer readings to the nearest 0.0001 fl. Tap mercury thermometers (for instance with a pencil) just before reading to avoid errors caused hy mercury sticking to the walls of the capillary. Ignite the charge at the start of the sixth minute and record the time, r, and the thermometer reading, r, (Note 6). Take the next two readings 0.5 min and 1 min after firing because of the rapid rise. Record subsequent readings at 1-min intervals on the minute until the difference hetween successive readings has been constant for 5 min. If the final temperature is above the bath temperature, the temperature rises to a maximum and then begins to fall. Record the time. f. and the thermumeter reading. In after the rate of change has become uniform.

Note δ —The calorimeter water temperature must be at the same temperature ($\pm 0.05^{\circ}$ C) for every determination, at the time of fining, if the Bureau of Mines method for radiation correction is to be used. See Annes A1.5.

10.3.7 Open the cover and remove the homb. Release the pressure at a uniform rate, such that the operation will require not less than 1 min. Open the homb and examine the humb interior. Discard the test if unburned sample or sooty deposits are found. Wash the interior of the bomb with distilled water containing the titration inducator until the washings are free of acid, and titrate the washings with standard sodium carbonate solution.

10.3.8 Remove and measure, or weigh, the combined pieces of unburned ignition (firing) wire and subtract from the original length or weight to determine the wire consumed in firing. If the wire is weighed, remove the ball of oxidized metal from the end of each piece of wire before weighing.

10.4 Cakulations:

10.4.1 Temperature Rise-Using data obtained as prescribed in 10.3.6, compute the corrected temperature rise, *s*, as follows:

1=1,-1;+6;+6;+6;

444

1 522000000214

*

AR301394

(1)

where:

1

- corrected temperature rise, "C or "F.
- = initial temperature reading at time of firł, ing.

1

- final temperature reading. ł.
- c: = thermometer, emergent stem correction. if required (see Note 7 and Annes A1.4).
- () = thermometer setting correction. if reguired (see Note 7 and Annex A1.3), and
- C. = radiation correction (see Note 7 and Annex A1.53.

Neny 7-With all mercury-in-glass thermometers. it is necessary to make corrections if the total calorific value is affered by 5 Bits or more. This represents a change of 0.0001°C or 0.002°F in a calorimeter using approximately 2000 g of water. Berkmann thermometen also require a setting correction and an emergent stent currection (Annes 41.3 and A1.4). Solid-stem ASTM Thermometers SAC and SAF do not require emergent stem corrections if all tests, including standardization, are performed within the same \$.5°C elu*F1 interval. If operating temperatures range beyond this limit, a differential emergent stem correction (Annev A1.4) must be applied to the corrected temperature mic. r. in all tests. including standardization. A radiation correction must be made in all tests, including standardization.

10.4.2 Thermachemical Corrections (see Appendix X1)-Compute the following for each tert:

- correction for the heat of formation of 6. 8 HNO., in Blu. Each millilitre of standard Na₂CO₂ is equivalent to 10.0 Biu. and
- c; = correction for heat of combustion of firing wire, in Blu (Note 8).
 - = 0.41 Biu/mm or 2.6 Biu/mg for No. 34 B & S gage Chromel C.
 - = 0.4% iu/mm or 3.2 Biu/mg for No. 34 B & S gage iron wire.

Nini K-There is no correction for platinum wire prevented the ignition energy is constant.

10.4.3 Compute the calorimeter energy equivalent, E, by substituting in the following:

$k = kh c + c_1 + c_1 V_1$

where:

Ŷ

7

• • •

E = calorimeter energy equivalent (Note 9),

- 11 = heat of combustion of heneoic acid, as stated in the National Bureau of Standards certificate. Blu/lb in gir,
- mass (weight in air) of benzoic acid, g.
- $r_1 = \text{titration correction (10.4.2)}$
- $r_1 = fusc wire correction (10.4.2), and$
- corrected temperature rise (10.4.1)

Nors 9-Using the units and corrections as given in 10.4.1 and 10.4.2. the energy equivalent of the colorimeter is such that the calonfic value of the coult sample will be obtained directly in British thermal write per pound when the mass sample is taken in grams. The units of the energy equivalent are therefore: (Biu/ Th) times (g/deg).

10.5 Repeat the procedure for a total of ten determinations. Compute the standard deviation as illustrated in Table 2.

11. Restandardization --

0 3288

11.1 Make checks on the energy equivalent value after changing the oxygen supply, after changing any part of the calorimeter, and at least once a month otherwise.

11.1.1 If a single new determination differs from the old standard value by 6 Biu/C (4 Biu/ F) the old standard value is suspect, thereby requiring a second test.

11.1.2 The difference between the two new determinations must not exceed \$ Btu/C (\$ Btu/ "H, and the average of the two new determingtions must not differ from the old standard by more than 4 Blu/C (3 Blu/F). If these requirements are met, do not change the calorimeter standard.

11.1 3. If the equirements given in 11.1.2 are not met, two more determinations must be run. The range of the four values must not exceed 14 Biu/C (8 Biu/F), and the average of the four new determinations must not differ from the old standard value by more than 3 Blu/C (2 Blu/ *F). If these requirements are met. do not change the calorimeter standard.

11.1.4 If the requirements given in 11.1.3 are not met. a fifth and sixth determination must be run. The range of the six new determinations must not exceed 17 Btu/C (10 Btu/F), and the average of the six values must not differ from the old standard value by more than 2 Btu/C (2 Btu/F). If these requirements are met. do not change the calorimeter standard.

11.1.5.10 the requirements given in 11.1.4 are not met, four more determinations must be run to complete a series of ten runs. The range of these ten results must not exceed 20 Btu/C 112 Btu/F), and the average of the ten new values must not differ from the old standard by more than 1 Blu/C (1 Blu/FL If these requirements are met, do not change the calorimeter standard,

11.1.6. If the requirements given in 11.1.5 are not met, the average value from the ten new

AR301395

(2)

445

values must be used for the new standard energy equivalent, provided that the standard deviation of the series does not exceed 6.5 Bug/C(3.6 Bug/ 4F).

11.2 The summary of the numerical requiremems at each stage of restandardization is given in Table 3.

12. Procedure for Coal and Coke Samples (Note 101

12.1 Thoroughly mix the analysis sample of coal or coke in the sample bottle and carefully weigh approximately 1 g of it into the sample holder. The sample shall be weighed to the nearest 0.0001 g. Make each determination in accordance with the procedure described in 10.3.2 through 10.3.8.

Nory 10-For anthracite, coke, and coal of high ash content, that do not readily burn completely, one of the following procedures are recommended (1). The inside of the sample holder is lined completely with ignited ashestos in a thin layer pressed well down into the angles, and the sample is then sprinkled evenly over the surface of the advestos. (2) The mass of the sample may be varied to obtain good ignition. If the mass is varied, it will be necessary to recalibrate the calorimeter so that the water equivalent will be based on the same temperature rise as that obtained with the new sample weight, (A) A known amount of beniew acid may be mixed with the sample. Proper allowance must be made for the heat of combustion of benzine acid when determening the caloritic value of the sample

Note 11-For the calorific value of coke, it is necwork to use 3 MPs (M atm) pressure for both stanto drawn and analysis.

12.2 Determine the sulfur content of the sample by any of the procedures described in Test Methods D 3177.

13. Calculations (Note 2)

 Compute the corrected temperature rise. 7, as shown in 10.4.1.

13.2 Thermachermical Convertions (Appendie ND—Compute the following for each test:

- $c_1 = correction$ for the heat of formation of HNOs, Bur Each milling of standard Błu.
- (c) = correction for heat of combustion of ignition wire. Blu.
 - = 0.41 Blu/mm or 2.6 Blu/mg for No. 14 B & S Chromel C wire.
 - = 0.49 Btu/mm or 3.2 Btu/mg for No. 34 B & S gage iron wire (Note 8), and

e. = correction for difference between heat of

formation of H₂SO₄ from the heat of formation HNO₁, Btu.

= 23.7 times percent of sulfur in sample tinies weight of sample, g

14. Cutorifie Value (Note 12)

14.1 Gross Calorific Value-Calculate the gross calorific value (gross heat of combustion at constant volume) Q. (gross) as follows:

$$Q_1(grow) = [(th) - c_1 - c_2 - c_1]/x$$
 (1)

where:

Ł

- Q. (gross) . gross calorific value Btu/lb. corrected temperature rise as calculated in 13.1.
 - = energy equivalent calculated in 10.4.3.

corrections as presented in 13.2. C1. C2. C3 and

mass of sample, g.

14.2 Net Calorific Value-Calculate the net coloritie value (net heat of combustion at a constant pressure). On (net) as follows:

$$O_{2}$$
 (net) = O_{1} (gross) = 10.30 (// x *)

where:

Q. (net) = net calorific value. Btu/lb.

Q. (gross) = gross calorific value, Buu/Ib, and 11 = total hydrogen, %.

Note 12-This calculation gives calonife value in Bus/lb: To obtain calorific value in J/g, see Appendix X ?.

15. Report

13.1 The results of the calorific value may be reported on any of a number of bases, differing from each other in the manner in mulsture is treated.

15.2 Use the percentage of monsture in the sample passing a No. 60 (250-µm) view. Test Method D 3173, to calculate the results of the analysis sample to a dry basis.

15.3 Procedures for converting the value obtained on the analysis sample to other bases are described in Method D 3180.

18. Prevision and Blue

16.1 The following criteria should be used for judging the acceptability of results 1983, probability) on split 60-mesh (250-jun) samples.

18.1.1 Repetitubility-Duplicate results by the same laboratory using the same operator and equipment, should not be considered suspect

446

 \int

unless they differ by more than 50 Btu/th on a, dry basis.

16.1.2 Reproducibility—The results submitted by two or more laboratories (different equipment, operators, date of test, and different portions of the same pulp) should not be considered

ر م

 suspect unless the two results differ by more than-100 Btu/lb on a dry basis.

16.1.3 Bias-There should be no bias because the equipment is standardized with a compound having a known heat of combustion.

1 ARI F 1 (l'aurille Vatur -
1 Big + Lissen J	4 Biu/th = 2 326 J/g
I Calinie a d 1865 J	1 # Big.16 # 1.0 cat/g

* International tables values:

TABLE 2 Sunders Deviations for Calerimeter Nanderdization*

Standardization	Culuma A	Colyma B	Column C
Numiter	Energy Equivalent (Miu/Rhi # ag/C)	(Code to 4400 (Column A-4400)	Column Br
	4412	12	144
j	4407	1	41
· j	4415	- 15	225
á l	4405	8	44
ġ	4404	▲ 1 1 1	16
i	4406		36
7	4404	. 4	61 · · · ·
ſ	4410	10	100
	4412	12	144
10 Stikt	<u>4404</u>	a.	420 420
vrags + 1 + 23, 10 + 192/101 +	43(E) = 43(W)		
maller + 1 - 1 - 1 - 1	26 idumn 1877/n) - 4dis - frag	r. tut = 10.4	
mdard destation • s • s distan	w = 10 d = 1 22		

" In this example the values of energy equivalent are typical for a calonmeter calibrated as thus, if the energy equivalent is multiplied by the temperature me in degrees Celous per gram of sample, the estimile value of the sample will be obtained in Bertich Thermal units per pound

TABLE 3 Kemmary of Numerless Regultements

Nint - Fest values exceeding table limits require additional

Number of Rons	Maumum Range of Koulta		Maximum Diffre- ence between X ₁ and X ₂ ^d	
	Bu/C	Biu/F	Ris/C	Bu/F
1		• • •	20	24
2		3	84	8)
. 4	14	8		22
6	12	10	#2	±2
ю	20	12		11

Values on this table have been musicled after stationed solvalation, and are therefore not precisely of a ratio from 2 B b + 4 D.

 ${}^{F}V_{s}$ — everage of original standard, and V_{T} = everage of clock runs

AR301397

ANNEX

(Mandatory Information) --

AL THERMOMETRIC CORRECTIONS-

A1.1 Thermometer Currections----

A1.1.1 It is necessary to make the following individual concerning if not making the corrections, would result in an equivalent change of 5.0 Big or more.

A1.1.2 Culibration Correction, shall be made in zesurdance with the calibration certificate furnished by the calibration authority.

ALI.3 Setting Correction. necessary for the Beckmann thermometer. It shall be made in accordance with the directions furnished by the calibration authority.

A1.1.4 Differential Emergent Stem Correction, the exlcutation of differential stem correction depends upon the way the thermometer was calibrated and how is is used. Two conditions are possible:

A1.1.4.1 Thermometers Calibrated in Total Immer-sion and Used in Partial Immersion-This emergent stem correction is made as follows:

Controlog =
$$C_{i} = K(\mu + i)(\mu + \mu - L - T)$$

where:

;÷. .

. ...

.,

4.

. .

- C. . emergent stem correction.
- 0.00016 for thermometers calibrated in "C. and •
- 0.00009 for thermometers calibrated in "F. scale reading to which the thermometer was im-
- mend.

mean temperature of emergent stem.

- initial temperature reading, and 1
- · final temperature reading. h

Nots At 1-Example Assume the point In to which the thermometer was immersed was 15°C; its initial reading, in was 24.127°C, its final reading, in was 27.876, the mean temperature of the emergent stem, T. was 26"C:

Differential stem correction, C, Differential stem correction, C, Diffusion (28 - 24) (28 + 24 - 16 - 26)

= 0.00(4°C).

A1.1.4.2 Thermoneters Calibrated and Used in Partial Immersion. Not at a Different Temperature than the Calibrated Temperature-This emergent sem correction is made as follows:

Conclusion =
$$C_{i} = X (i_{i} = i_{i}) (i_{i} = i_{i})$$

where.

C. . emergent stem correction.

- 0.00016 for thermometers calibrated in "C. and
 - U.IXXIV for thermometers calibrated in "F. · initial temperature reading.
- L · final temperature reading.

30

- · observed stem temperature, and
- sizm temperature at which the thermometer was calibrated.

Note A1.2-Example: Assume the initial reading. r, was 80°F, the final reading, r. was \$5°F, and that the observed stem temperature, to was \$2°F, and the calibration temperature, to was 72°F, then

Differential stem correction. C. - 0.00009 (16 - 10) (82 - 72)

= 0.005°F.

A1.1.5 Radiation Correction-These are used to culculate heat loss to the water jacket. They are based on the Dickinson formula," the Regnault-Plaundler formula," or the U.S. Bureau of Mines method." The same method of determining the radiation correction must be used consistently in calibration and test measurements.

A1.1.5.1 Dickinson Formula:

where:

۶.

1

C.

Ħ

S

C. = radiation correction.

- rate of rise in temperature per minute in the preliminary period. rate of rise of temperature per minute in the final
 - period (if temperature is fulling, of is negative). .
 - firing temperature. .
 - final temperature, being the first temperature after which the rate of change is constant.
 - time at temperature, I, min.
- h . time at temperature. $t_i \neq 0.50$ ($t_i = t_i$). min. and
 - = sime at temperature. In min. A1.1.5.2 Regnault-Phaumiller Formula
 - $C_i = m_i + kS$

where:

- rediation correction. . number of minutes in the com-
- busiton period,
 - = (1 = 1)/(1* = 1°).
 - = in 1 + (1/2) (i, + in ne".
 - = sverage temperature during the preliminary period.
- average temperature during the final period.

1. 1. 1.

*Dickensin, H. C., Bulletin U.S. Bureau of Mandards, Not 11, 1451, p. 189

Prophyles, L., January der Physick (Leiperg), ANPA And

129 (John p. 102). **Methods of Analyzing and Testing Cost and Color* E.S. Amount Methods Bulliour 8(8), XMM14, 196*, pp. 16-1*

44%

during the combustion period. and

1. - 1

= sum of $t_1, t_2, t_3, \dots, t_n = 1$.

A1.1.5.3 Bureau of Mines Method-A table of m. diation correction can be established so that only the initial and final readings are required to determine the heat value of any fuel. This may be done by carrying out a series of tests utilizing the procedure described in Section 10, using the following conditions. Regulate the annunt of sample burned so that a senes of determinations is made in which different temperature rises are obtained. For all determinations, keep the water

jacket temperature constant. fire the bomb at the same initial temperature, and have the same time. f + 1. clapse 122 st between the initial and final readings Determine the radiation corrections for each of the series of temperature rises using the Dickinson method (see A1 1.5.1), or the Regnault-Plaundler method (see A1.1.5.2). These corrections are constant for a given temperature rise. From the series of readings a table or sraph is plotted to show radiation correction versus temperature rise. Once the table of graph is established. the radiation corrections can be obtained from it until there is a major change in the equipment.

APPENDIXES

(Nonmandatory Information)

XI. THERMOCHEMICAL CORRECTIONS-

X1.1 Energy of Formation of Nitric Acid—A correction. e. (10.4.2 and 13.2) is applied for the acid inration. This correction is based on the assumptions (1) that all the acid titrated is HNO, formed by the following reaction: 12 Nr (g) + % Or (g) + % HrO(1) = HNOr (in SOU mot H2O), and (2) that the energy of formation of HNO. in approximately 500 mol of water under bornh runditions is minus 39.0 LJ/mol.16

X1.1.1 A convenient concentration of Na;CO, 18 0.394 N (20.9 g Na; CO_{3} /1000 m1.3 which gives $c_{1} = 10$ multiply 1', where 1' is the volume of Na; CO_{3} in millilitres. The factor 10.0 (0.344 x 59.0 + 2.326) is to be used for esteulating calorific value in Blu/lb. For other units see Table X2.1. When H₂SO₆ is also present. a part of the correction for HISO, is contained in the er correction and the remainder in the er correction.

X1.2 Energy of Formation of Sulfurn: Acid-By definition (ace Definitions D 121) the gross calorific value is obtained when the product of the combustion of sulfur in the sample is SO2 (g). However, in actual tomb combusion processes, all the sulfur is found as HisO, in the bomb washings. A correction of thee 13.21 in applied for the sulfur that is converted to HisOs. This correction is based upon the energy of formation of H₂SO, in solutions, such as will be present in the homb at the end of a combustion. This energy is taken as -295.0 W/mol." A correction of 2 times 54.0 W/mot of sulfur was applied in the ri correction, so the additional correction necessary is 295.0 = (2 times \$9.0) = 177 LJ/mol, or \$ \$2 LJ/g of sulfur in the sample (\$5.2 I times weight of sample in grams times percent sulfur in sample). This causes of to be 23,7 times weight of samples in grams times percent sulfur in sample. The factor 23.7 (= 55.2/2.326) for e, (see 13.2) is to be used for calculating exterilic value in Biu/th. For other units. see Appendix X2. The values showe are based on a coult containing shout 5% sulfur and shout 5% hydrogen.

ويوازيلون بالجوت وكرافته فأبرز

The assumption is also made that the HSO, is dissolved entirely in the water condensed during combustion of the sample.¹²

X1.2.1 If a 1-s sample of such a fuel is byrned, the resulting H2SO2 condensed with water formed on the walls of the bomb will have a ratio of about 15 mill of water to 1 mot of HisO. For this concentration the energy of the reaction SO2 (g) + 12 O2 (g) + H10 (1) = 11:50. (in 15 mol of H₁O) under the conditions of the homh process is -295 LJ/mol." Basing the calculation upon a sample of comparatively large sulfur content reduces the possible overall errors, because, for small percentages of sulfur, the correction is smaller,

X1.3 Fust (leninm) Wire-Calculate the energy contributed by burning the fuse wire in accordance with the directions furnished by the supplier of the wire For example, the energy of the combustion of No. 34 B & S gape Chromel C wire is 6.0 J/mg or approximatch UNS J/mm. For calculating ra for use in Eqs 2 and 3. these give rg = 0.41 times length (mm) of wire or 1y = 2.6 times weight (mg) of wire. The energy required to melt a platinum wire to constant for each experiment if the same amount of platinum wire is used. As the energy is small, its effect is essentially cancelled out in the relationship between the standardwation experiments and the calorific value determinations, and it can be neglected. The factors listed above far of \$10.4.2 and 1.3.2) are suitable for calculating culurific value in Bus/Ib. For other units, see Appendix X2.

⁴ Calculated from data in National Noreau of Mandarch Sochrings Note 27(1) 5 ¹¹ Calculated from data in National Burgar of Mandarch

Concular NB: ¹⁴ Marti B. A., and Parker C., "Nucleos on Bornh Calorimeter UNE formation of Sulfary, Acid." J and J CLEB, Vol. 17, 1949. 8 171

X2. REPORTING RESULTS IN OTHER UNITS

1

X21 Reporting Results in Juniles per Gram

X2.11 The gross calorifie value can be expressed in joules per gram, calones per gram, or British thermal units per pound. The relationships between these units are given in Table I. X2.1.2 Because the energy of combustion of the

reference material is measured and certified by the National Bureau of Standards in joules per gram, the most straightforward usage of the reference material would lead to the calorific value of the fuel in joules per yrum. To camy out this procedure, we make changes outlined in N2.1.3 through X2.1.5.

X2 1 * For calculating energy equivalent, substitute= es', ep', and es' = corrections, with units of joules (see Eq.2: 12NP X2 1).

$$E = [I (II'_{2}) + r_{1}' + r_{2}']/t \qquad (2')$$

= heat of combustion of reference material, with units of joules per

gram weight in air 13/g from the

certificate for the NBS bensore

where the meanings of the symbols in Eu 21 are the sinclas in Eq.2 except that:

acid) and

- energy equivalent with units of joules per temperature unit.

Ľ 11

Antoine and a second second

2

ы.-А.--

-..... o', rj', and rj' = corrections, with units of joules (see-Table X2.11

X2.1.4 For calculating gross calorific value, substitute Eq 3' for Eq 3'

$$Q_{1}(gross) = \{(tE') - e_{1}' - e_{2}' - e_{3}'\}/g = \{(N)\}$$

where the meanings of the symbols in Eq.3" are the same as in Eq 3 except that:

Q. (genu) · gross culorific value with units of joules per gram (weight in air).

= energy equivalent, with units of

AR301400

joules per temperature unit, and

same laboratory, using the same opirator and equipment, should not be considered suspect unless they differ hy more than 120 J/g

X2.1.5.2 Reproducibility-The results submitted by two or more fuboratories tdifferent equipment, operators, date of test, and different portions of the same pulp) should not be considered suspect unless the results differ by more than 240 J/g.

TABLE N2.1 Abstrative Thermachemical Correction-Factors (L'aits in Jouins)⁴

Correction	Multiplication Factor	Multiply By
64° 111NO41 65° 13458041	20 J/mL 55.2 J/cgS	mL of 0.34 N NagCOs proventage of solitor in sample tumes mass of sample in promo
iy' (lure wite)	194 1/ 1/10	Irright constant No. 34 B & S gap Chaund C way
in 17' Hune week 1	1.14 J/mm	length terms of No. 34 B & S
ay" (lung wate) af	34J/mg	main (ma) of Chromet C wire
iy' thue water	.7.4 3/mg	musi ting) of iron wire

* for he used in Fus 2" and 3" only

The some on Source in Ferring and Materials cakes no position requesting the valuates of any polent rights asserted in connection with any new memory on this sepadard. Users of this seadard are expressly advised that desemblation of the valuary of our such parent rights, and the risk of inforgement of such rights, are entirely them our regionsolvery.

This sumdard is subject to resisting at any time by the responsible technical committies and must be restered every the secon and if the second either trapponed to although a Your comments are instead either for second at this standard or two additional standards and should be addressed to ASTM Headquarter. Your comments will reserve careful consideration at a moving of the symmittle technical committees which is a mass attend. If you had that your comments will reserve an overed a two heaving of the responsible technical committees which is a mass attend. If you had that your comments have me recovered a two heaving on chald make some second most to the ASTM Committee on Standards, 1916 Rule 50. Philadelphia, PA 19103.

gO. 15 g of powdered K2SO4. f of H₂SO₄. Place the flask in an ning and heat gradually. Then boil t solution clears. Continue boili...ional 30 min. Cool. add about ater. cool below 25°C; add 25 mL ting and mix to precipitate the d pinch of zinc to prevent bumpfinsk. and add a layer of NaOH fficient to make the contents d e. Do not agitate the mixture. meet the flask to the digestion denser. Have the tip of the conrsed in the H₂BO₂ solution (this n! sured) in the receiver and then o mix the contents thoroughly. ie ammonia has distilled (at least stillate). Titrate with standard acid irl icator to violet end point.

ate the percentage nitrogen as fol-

 $\dots ed), \% = (A \times B) \times 0.14$

rel of 0.1 N to 0.3 N H2SO4 retitration of the solution, and lity of the H2SO.

the nearest 0.1% the nitrogen : as-received sample.

i Accuracy

ision of this method has not aed. Data are being sought that le "or use in developing precision

be validity of any patent rights asserted in sty advised that determination of the validity • l • responsibility.

where and must be reviewed every five years Where are mass or revenue or for additional or revision of this standard or for additional the careful consideration at a meeting of the e not received a fair hearing you should Aug. Pa. 19103.

Designation: D 2974 - 84

Vol. 04.08

Standard Test Methods for MOISTURE, ASH, AND ORGANIC MATTER OF PEAT MATERIALS'

This standard is issued under the fixed designation D 2974: the number immediately following the designation indicates the year of original adoption or. in the case of revision, the year of last revision. A number in parentheses indicates the year of last reapproval A superscript epsilon (+) indicates an editorial change since the last revision or reapproval.

1. Scope

1.1 These test methods² cover measurement of the weight percentage of moisture, ash, and organic matter in peat materials, including moss, humus, and reed-sedge types.

1.2 This standard may involve hazardous materials, operations, and equipment. This standard does not purport to address all of the safety problems associated with its use. It is the responsibility of whoever uses this standard to consult and establish appropriate safety and health practices and determine the applicability of regulatory limitations prior to use.

2. Summary of Methods -

2.1 Method I-Moisture is determined by drying a peat sample at 105°C. The loss in weight represents the total moisture content expressed as percentage by weight of the as-received sample.

2.2 Method II-This is an alternative moisture method which removes the total moisture in two steps: (1) evaporation of moisture by air at room temperature (air-drying) and (2) the subsequent loss in weight of the air-dried sample at 105°C. This method provides a more stable sample, the air-dried sample, when tests for nitrogen, pH, cation exchange, sand content, etc., are to be made.

2.3 Ash content of a peat sample is determined by igniting the oven-dried sample from the total moisture determination in a muffle furnance at 550°C. The substance remaining after ignition is the ash, and this includes mineral impurities such as sand. The ash content is expressed as a weight percentage of the oven-dried sample.

2.4 Organic matter is determined by difference.

3. Apparatus

3.1 Oven, regulated to a constant temperature of 105 ± 5°C.

3.2 Muffle Furnace, regulated to a constant temperature of 550°C.

3.3 Evaporating Dishes. high silica or porcelain. not less than 75-mL capacity.

3.4 Blender, high-speed.

3.5 Aluminum Foil, heavy-duty.

3.6 Porcelain pan, spoons, etc.

4. Preparation of Sample

4.1 Place a representative field sample on a square rubber sheet, paper, or oil cloth. Reduce the sample to the quantity required by quartering and place in a moistureproof container. Work rapidly to prevent moisture losses.

MOISTURE

5. Procedure-Method I

5.1 Mix the sample thoroughly and place a test specimen of at least 50 g in an ignited and weighed (with fitted heavy-duty aluminum foil cover) high silica or porcelain evaporating dish of not less than 100 mL capacity. Crush soft lumps with a spoon or spatula. Cover immediately with the aluminum foil and weigh to the nearest 0.01 g. Dry uncovered for 16 h at 105°C. Remove from the oven, cover tightly, cool. and weigh.

These test methods are under the junsdiction of ASTM Committee D-18 on Soil and Rock and are the direct responsebility of Subcommittee D18.18 on Peats and Related Materials.

Current edition approved Oct. 26, 1934. Published December 1984. Originally published as D 2974 - 71. Last previous edition D2974-71. These test methods A Brash and the and an extensive

review by Committee D-18. -



where:

6. Calculation~

system is being used.

7. Procedure-Method II

where:

minations.

lows:

Equivalent sample weight.

and proceed as in Method I.

8. Calculation

.

6.1 Calculate the moisture content as follows: Moisture: $\% = [(A - B) \times 100]/.4$

NOTE-In geotechnical engineering the moisture

7.1 Use Method II when pH. nitrogen, sand

content, cation exchange capacity, etc., are to be

tested. Mix thoroughly and weigh a 100 to 300g representative sample and spread evenly on a

large flat pan. Crush soft lumps with a spoon or

spatula and let the sample come to moisture

equilibrium with room air, not less than 24 h.

Stir occasionally to maintain maximum air ex-

posure of the entire sample. When the weight is constant, calculate the loss in weight as percent

moisture removed by air drying. Grind a representative portion of the air-dried sample 1 to 2

min in a high-speed blender. Use the ground

portion for moisture, 2sh, nitrogen, etc., deter-

7.2 Thoroughly mix the air-dried, ground sample and weigh to the nearest 0.01g the equiv-

alent of 50 g of test specimen on the as-received basis. Determine the grams of air-dried sample

equivalent to 50 g of as-received sample as fol-

Place the weighed sample in an ignited and

weighed (with fitted heavy-duty aluminum foil cover) high silica or porcelain evaporating dish

8.1 Calculate the moisture content as follows:

g = 50.0 - [(50 × percent moisture removed

content is commonly expressed as a proportion of the oven-dried mass. Care must be taken to indicate which

A =grams of as-received test specimen, and

B = grams of oven-dried test specimen...

B =grams of oven-dried sample.

ASH

Moisture, $\% = (50 - B) \times 2$

9. Procedure-

9.1 Place an uncovered (retain cover for weighing) high silica or porcelain dish containing the dried test specimen from the moisture determinations in a muffle furnace. Gradually bring to 550°C and hold until completely ashed. Cover with the retained aluminum foil cover, cool, and weigh.

10. Calculation

10.1 Calculate the ash content as follows:

Ash. % = (C × 100)/B

where:

C = grams of ash. and

B = grams of oven-dried test specimen.

ORGANIC MATTER

11. Procedure

11.1 Determine the amount of organic matter by difference, as follows:

Organic matter, % = 100.0 - D

where:

D = weight percentage of ash.

REPORT

12. Report

12.1 Report all results to the nearest 0.1 %.

12.2 Indicate whether moisture contents are by proportion of as-received mass or oven-dried mass.

13. Precision and Bias

13.1 The precision and bias of these test methods have not been determined. Data are being sought for use in developing a precision and bias statement.

The American Sectors for Texting and Materials takes no position respecting the validity of any patent rights asserted in connection with any tiem mentioned in this standard. Users of this standard are expressly advised that determination of the validity of any such purent rights, and the risk of infringement of such rights, are entirely their own responsibility

in air drying)/100]

This standard is subject to revision at any time by the responsible technical committee and must be reviewed every five years and it we revised, either reapproved or willdrawn. Your comments are invited either for revision of this standard or for additional undurity and should be addressed to ASTM Headquarters. Your comments will receive careful consideration at a meeting of the responsible technical commutee which you may attend. If you feel that your comments have not received a fair hearing you should make whit views known to the ASTM Committee on Standards, 1916 Race St., Philadelphia, Po. 19103


498

Standard Test Meth SAND CONTENT

Last ASTM Designations

This test method covers determi peat.

Formerly under the jurisdiction This test method was withdrawn

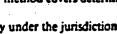

. CARBON, TOTAL ORGANIC AND INORGANIC

Carbon may exist in sediment and water samples as either inorganic or organic compounds. Inorganic carbon is present as carbonates, bicarbonates, and possibly free carbon dioxide. Specific types of compounds that are considered to be included in the organic carbon fraction are nonvolatile organic compounds (sugars), volatile organic compounds (mercaptans), partially volatile compounds (oils), and particulate carbonaceous materials (cellulose).^{1,2[#]}


The basis of the method is the catalytic or chemical oxidation of carbon in carbon-containing compounds to carbon dioxide followed by the quantification of the carbon dioxide produced. Alternately, the carbon may be reduced to methane and appropriately quantified. It follows, then, that the distinction between inorganic carbon and organic carbon is the method of sample pretreatment. There are presently two procedures for defining this separation. One method is based on sample treatment with a strong acid. Analysis of an untreated sample is a measure of total carbon while analysis of the acid-treated fraction is a measure of organic carbon. Inorganic carbon is calculated by subtraction. The second method of separation is based on differential thermal combustion with organic compounds being converted to carbon dioxide at 500° C to 650° C^{3,4} and inorganic carbon being converted to carbon dioxide at 950° C to 1300° C.^{4,5}

Sample Handling and Storage

: ÷.

Flowcharts for the handling of samples intended for organic carbon and inorganic carbon analysis are presented in Figure 3-6 and Figure 3-7. Water and sediment samples to be analyzed for inorganic carbon may be stored in glass or plastic containers. There is no effective preservative because of the carbon dioxide reserve in the atmosphere. The only precaution that can be taken for inorganic

* References for this procedure can be found on page 3-76. AR301403

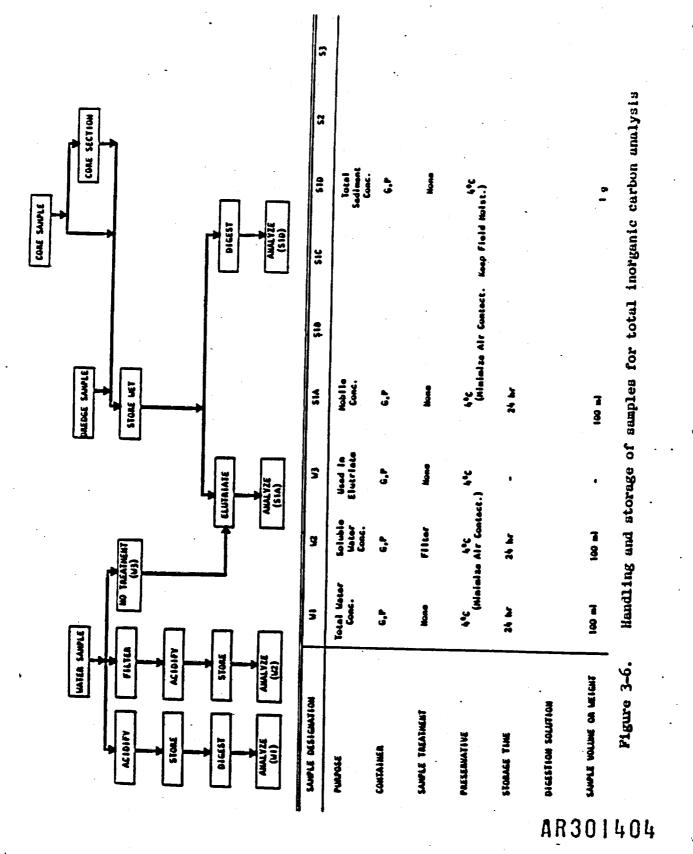

	\$10 \$2 \$3	Total Sediment	conc. conc. conc. G,P G,P G,P	None Alr dry Freeze	4°C None None	•	0.5-5.0 9 0.5-5.0 9 0.5-5.0	total organic carbon analysis
	AJ SIA	Bsed in Mobile Elutriate Conc.	£. 3	Rone	4.C 4.C	7 7 7	Var læle	Handling and storage of samples for total organic
	24	Total Water Soluble B3 Conc. Nater Elle Conc. Conc.	. U	Mane fliter A	MC1 or M ₃ 50, MC1 or M ₃ 50, 4 to p ^{M2} to p ^{M2} 4°C 4°C	24 hr	30-100 m) 30-100 ml Var	Handling and store
		1	CONTAINER	SAMPLE TREATMENT		AR30		Figure 3-7.

ŧ.

建設

carbon is to completely fill the sample container at the time of sampling (exclude all air bubbles), tightly seal the container, and complete the analysis immediately (Figure 3-6).

Water samples for organic carbon analysis should be stored in glass containers unless substitute containers have been shown not to affect total organic carbon (TOC) analyses. Samples should be processed as soon as possible (within 24 hr if possible) to minimize change due to chemical or biological oxidation. Atmospheric uptake of carbon dioxide is less critical since it would be evolved when the sample is acidified prior to analysis. Sediment samples for organic carbon analysis may be stored in either plastic or glass containers (Figure 3-7). Air drying of sediments (S2) may lead to low TOC results due to oxidation or volatilization. Therefore, moist storage (SID) or frozen storage (S3) would be the preferred method of storage. If samples are frozen, excessive temperatures should not be used to thaw the samples.

3-68

AR301406

Procedure for Water Samples (W1, W2, S1A)

Method 1: Infrared Analysis ",7

Apperatus

Sample homogenizer such as a Waring blender or ultrasonic blender

Magnetic stirrer

Hypodermic syringe

Total carbon analyzer, either a single channel or a dual channel instrument (Dow-Beckman Carbonaceous Analyzer Model No. 915, Dohrmann Envirotech DC-50 carbon analyzer, Oceanography International Total Carbon Analyzer, Leco, or equivalent)

Reagents

Distilled water: the distilled water used in the preparation of standards and dilution of samples should be of the highest quality in order to have a small blank.

Organic carbon, stock solution, 1000 mg/l C: dissolve 2.125 g anhydrous potassium biphthalate, KHC: H+O+, in distilled water and dilute to 1 ^l in a volumetric flask.

Organic carbon, standard solutions: prepare standard solution by dilution of the stock solution as required.

Inorganic carbon, stock solution, 1000 mg/2: dissolve 3.500 g sodium bicarbonate, NaHCO3, and 4.418 g sodium carbonate, Na2CO3, in distilled water in a 1-2 volumetric flask and make up to the mark.

Inorganic carbon, standard solution: prepare standards from the stock solution as required.

Packing for total carbon tube: dissolve 20 g cobalt nitrate, $CO(NO_2)_2$ · GH_2O_1 in 50 ml distilled water. Add this solution to 15 g long-fiber asbestos in a porcelain evaporating dish. Mix and evaporate to dryness on a steam bath. Place the dish in a muffle furnace and bring to 950°C. After 1 to 2 hr at this temperature, remove the dish and allow to cool. Break up any large lumps and mix adequately but not excessively. With the combustion tube held in a vertical position, taper joint up, put about 1/2 in. of untreated asbestos in the tube first, then transfer in small amounts, approximately 1 g of catalyst into the tube with forceps or tweezers. As it is added, tap or push the material gently with a 1/4-in. glass rod. Do not force the packing. The weight of the rod itself is sufficient to compress the material. When completed, the length of the packing should be about 6 37 | 407 6 cm. Test the packed tube by measuring the flow rate 37 gas through it at room temperature, and then at 750°C. The rate should not drop more than 20 percent.

Packing for carbonate tube (dual channel instrument): place a small wad of quartz wool or asbestos near the exit end of the carbonate evolution tube. From the entrance end add 6 to 12 mesh quartz chips, allowing these to collect against the wad to a length of 10 cm. Pour an excess of 85 percent phosphoric acid, H₃PO₄, into the tube while holding it vertically and allow the excess to drain out.

Nitrogen gas, carbon dioxide free.

Procedure

Turn on the infrared analyzer, recorder, and tube furnaces, setting the total carbon furnace at 950°C and the carbonate furnace at 175°C. Allow sufficient warm-up time for stable, drift-free operation; about 2 hr is required. If used daily, the analyzer can be left on continuously. Adjust the oxygen flow rate to 80 to 100 ml/min through the total carbon tube. With other instruments, follow manufacturer's directions to warm up the instrument.

Immediately prior to carrying out calibrations or analyses, inject several portions of the appropriate standard into the tube to be used, until constant readings are obtained. The actual injection technique is as follows: rinse the syringe several times with the solution to be analyzed, fill, and adjust the volume to be pipeted. Wipe off the excess with soft paper tissue, taking care that no lint adheres to the needle. Remove the plug from the syringe holder, insert the sample syringe, and inject the sample into the combustion tube with a single, rapid movement of the thumb. Leave the syringe in the holder until the flow rate returns to normal, then replace it with the plug.

Successively introduce a convenient sized aliquot (20 to 50 µl) of each organic carbon standard and a blank into the total carbon tube and record peak heights. Between injections allow the recorder pen to return to its baseline. When a dual channel instrument is used, the standardization procedure must be repeated using carbonate standards to calibrate the low temperature channel.

Thoroughly mix the sample. Inject a convenient sized aliquot (20 to 50 μ 1) of the sample into the total carbon tube and AR301408

record the peak height. This result is a measure of the organic carbon concentration and the inorganic carbon concentration of the sample.

Thoroughly mix the sample using a Waring blender or an ultrasonic homogenizer. Transfer 10 to 15 ml of sample to a 30-ml beaker and acidify with concentrated HCl to a pH of 2 or less. Furge the sample with carbon dioxide free nitrogen gas for 5 to 10 min. Plastic tubing should not be used during the purging process unless it has been previously shown that it will not add organic carbon to the sample.

Mix the acidified sample on a magnetic stirrer. While stirring, withdraw a subsample from the beaker using a hypodermic needle with a 150-µm opening. Inject the sample into the carbon analyzer to be used and record the peak height. This result is a measure of the organic carbon concentration of the sample.

Using either clear or filtered water samples, analytical precision will approach 1 to 2 percent or 1 to 2 mg/l carbon, whichever is greater. Analytical precision for unfiltered water samples will increase to 5 to 10 percent because of the difficulty associated with sampling particulate matter and the fact that the needle opening of the syringe limits the maximum size of the particles that can be included in the sample.

Calculations

. .

-

. .

-

<u>Dual-channel instrument.</u> Prepare calibration curves derived from the peak heights obtained with the standard total carbon and inorganic carbon solutions.

Determine the concentration of total carbon and inorganic carbon in the sample by comparing sample peak heights with the calibration curves.

Determine the concentration of total inorganic carbon in the sample by subtracting the organic carbon value from the total carbon value.

Single-channel instrument. Prepare a calibration curve derived from the peak heights obtained with the standard total carbon solutions. Determine the total carbon concentration in the sample by comparing the peak height of the first sample injection with the AR301409 calibration curve. Determine the organic carbon concentration in the sample by comparing the peak height of the second sample injection with the calibration curve. Inorganic carbon concentrations are calculated by subtracting the organic carbon concentration from the total carbon concentration.

Procedures for Sediment Samples (S1D, S3)

Method 1: Sample Ignition

Apparatus

Induction furnace such as the Leco WR-12, Dohrmann DC-50, Coleman CH analyzer, or Perkin Elmer 240 elemental analyzer

Combustion boats

Microbalance

Desiccator

Reagents

10 percent hydrochloric acid: mix 100 ml concentrated HCl with 900 ml distilled water.

Copper oxide fines.

Benzoic acid.

Procedure

Dry at 70°C and grind the sediment sample.

Weigh a combustion boat and record the weight. Place 0.2 to 0.5 g homogenized sediment in the combustion boat and reweigh. Combustion boats should not be handled with the bare hand during this process.

If total carbon or inorganic carbon is to be determined, Cupric oxide fines may be added to the sample to assist in combustion. Combust the sample in an induction furnace. Record the result as total carbon.

If organic carbon is to be determined, treat a known weight of dried sediment with several drops of 10 percent HCL. Wait until the effervescing is completed and add more acid. Continue this process until the incremental addition of acid causes no further effervescence. Do not add too much acid at one time as this may cause loss of sample due to frothing.

Dry the sample at 70° C and place in a desiccator. Add Cupric oxide fines, combust the sample in an induction furnace, and record the result as organic carbon.

3-73

Calculations

The carbon content of the sample can be calculated as:

SC = Weight of tube (after-before) x 27.29 sample weight

Derivation of factor:

27.29 = 12.011 (molecular weight carbon) * 100%

When the total sample results are used, the result is percent carbon in the sample. When acid-treated samples are used, the result is percent organic carbon. Inorganic carbon is calculated as total carbon minus organic carbon.

Method 2: Differential Combustion ", 5

Apparatus

Sargent programmed microcombustion apparatus or equivalent Microbalance

Procedure

4...

÷.

Air dry the sediment sample. Using a mortar and pestle, grind the sample to pass a 100-mesh screen.

Combust a known weight of sediment at a programmed heating rate of 300° to 950° C in 10 min and then maintain 950° C for 20 min. Trap the CO₂ in ascarite and record the weight as total carbon. A sample size should be selected that will produce 25 to 50 mg CO₂.

Weigh a second portion of the dried sediment. Combust this sample at a programmed rate of 300° to 650° C in 10 min and maintain 650° C for 20 min. Trap the CO₂ in ascarite and record the weight as organic carbon.

Calculations

The total carbon concentration, C_t , of the sample (in mg/g) is calculated as follows:

$$C_t = \frac{\binom{x}{t} \binom{12}{44}}{\binom{g}{5}}$$

where

xt = weight of CO₂ evolved at 950°C, mg

3-74

g = weight of sample combusted, g

The organic carbon, C_0 , concentration of the sample (in mg/g) is calculated as follows:

$$C_{0} = \frac{(x_{0})(\frac{12}{62})}{(g)}$$

where

11

 x_{o} = weight of CO₂ evolved at 650°C, mg

g = weight of sample combusted, g

Inorganic carbon,
$$C_{I}$$
, (in mg/g) is calculated as:
 $C_{I} = C_{t} - C_{o}$

Method 3: Wet Combustion^{4,8}

A third method has been used for carbon in sediments. This is based on the oxidation of the sample with dichromate and back titration of the sample with ferrous ammonium sulfate. References are provided for the procedure but details are not given. The procedure is similar to the chemical oxygen demand test which is not specific for carbon. The wet combustion method is a redox procedure and any reduced chemicals in the sediment samples (ferrous iron, manganous manganese, sulfide) will react with the dichromate. Therefore, this procedure is not recommended unless other instrumentation is not available.

References

- 1. U. S. Environmental Protection Agency. "Manual of Methods for Chemical Analysis of Water and Wastes." Methods Development and Quality Assurance Research Laboratory, National Environmental Research Center; Cincinnati, Ohio. 298 p. (1974).
- 2. U. S. Environmental Protection Agency. "Methods for Chemical Analysis of Water and Wastes." Environmental Monitoring and Support Laboratory, Office of Research and Development, IPA; Cincinnati, Ohio (1979).
- 3. Giovannini, G., Poggio, G., and Sequi, P. "Use of an Automatic CHN Analyzer to Determine Organic and Inorganic Carbon in Soils." Unpublished Report, Laboratory of Soil Chemistry, via Corridoni, Pisa, Italy. 9 p. (1975).
- Konrad, J. G., Chesters, G., and Keeney, D. R. "Determination of Organic- and Carbonate-Carbon in Freshwater Lake Sediments by a Microcombustion Procedure." J. Thermal Analysis 2:199-208 (1970).
- Kemp, A. L. W. "Organic Matter in the Sediments of Lakes Ontario and Erie." Proc. 12th Conference Great Lakes Research 12:237-249 (1969).
- 6. Environment Canada. "Analytical Methods Manual," Inland Waters Directorate, Water Quality Branch; Ottawa, Canada (1974).
- 7. American Public Health Association. <u>Standard Methods for the</u> <u>Examination of Water and Wastewater</u>. APHA; New York, New York. 1193 p. (1976).
- Gaudette, H. E., Flight, W. R., Toner, L., and Folger, D. W.
 "An Inexpensive Titration Method for the Determination of Organic Carbon in Recent Sediments." J. Sed. Petrology 44:249-253 (1974).

3-76

553 09 04

AR301415

REM III SAS NO.

Ebasco Services Incorporated REM III Program Zone Program Management Office 2000 15th Street North Arlington, VA 22201 Phone: (703) 558-7555

> REM III TEAM LABORATORY SPECIAL ANALYTICAL SERVICES Regional Office Request

Regional	Transmittal

Telephone Request

A. Regional Office: <u>NUS-Pittsburgh</u> and Site Name: <u>Croydon TCE</u>

B. Regional Laboratory Service Coordinator: Greg Zimmerman

C. Telephone Number: (412) 788-1080

D. Date of Request: September 22, 1988

Please provide below a description of your request for Special Analytical Services under the REM III Program. In order to most efficiently obtain capability for your request, please address the following considerations, if applicable. Incomplete or erroneous information may result in delay in processing your request. Provide response on additional sheets or attach supplementary information as needed. 1. General description of analytical service requested: <u>Analysis of 7</u> sediment samples for TOC, and Grain Size Distribution.

 Definition and number of work units involved (specify whether whole samples or fractions; whether organics or inorganics; whether aqueous or soil and sediments; or whether low, medium, or high concentration: <u>7 low concentration sediment samples.</u>

3. Purpose of analysis (specify whether Remedial or Enforcement Action, RCRA, NPDES, etc.): <u>RI/FS</u>

4. Data Quality Objectives (specify whether screening analysis, engineering support or confirmational analysis): <u>Engineering Support</u> <u>DOO Level III.</u>

5. Estimated date(s) of collection: November 14 - 18, 1988

AR301416

- Estimated date(s) and method of shipment: <u>November 14 18, 1988.</u> <u>Samples will be shipped by overnight air carrier. These dates are</u> <u>tentative and are dependent on the project remaining on schedule.</u>
- Approximate number of days results are required after receipt of samples by laboratory: <u>40 days</u>
- Analytical protocol required (attach copy if protocol other than currently approved under this contract or CLP program):
 <u>Grain Size Distribution - ASTM D 422-63</u>
 <u>TOC - page 3-73 from "Procedures for Handling and Chemical Analysis of</u>
 <u>Sediment and Water Samples</u>", EPA, May 1981.

All methods are attached.

a.

- 9. Special technical instructions (if outside protocol requirements, specify compound names, detection limits, etc.):_____
- 10. Analytical results required (if known, specify format for data sheets, QA/QC reports, chain-of-custody, documentation, etc.). If not completed, format of results will be left to discretion of the laboratory: <u>Data sheets, chain-of-custody forms, all calculations,</u> and SAS Request Form.
- 11. Other instructions (use additional sheets or attach supplementary information, as needed):_____

12. Name of sampling/shipping contact: Greg Zimmerman

Phone: (412) 788-1080

ð.

, 4.-A **a** -

DATA REQUIREMENTS

I.

· · ·	mg/l	PRECISION DESIRED			
PARAMETER	DETECTION LIMIT	(+ % or Conc.)			
		en de la construcción de la constru La construcción de la construcción d			
Grain Size	Not Specified	Not Specified			
TOC	Not Specified	Not Specified			

II. OUALITY CONTROL REQUIREMENTS

SAMPLES	CONTROL	PRECISION DESIRED			
CONTROL REQUIRED	FREQUENCY OF SAMPLES				
		•			
None Specified					
	• • •	•			
) 					
· · · · · · · · · · · · · · · · · · ·					

III. ACTION REQUIRED IF LIMITS ARE EXCEEDED:

None Specified

M Designation: D 422 - 63 (Reapproved 1972)⁴¹

Standard Method for Particle-Size Analysis of Soils¹

This standard is issued under the fixed designation D 422; the number immediately following the designation indicates the year of original adoption or, in the case of revision, the year of last revision. A number in parentheses indicates the year of last reapproval. A superscript epsilon (+) indicates an editorial change since the last revision of reapproval.

"Note-Section 2 was added editorially and subsequent sections renumbered in July 1984

1. Scope-

1.1 This method covers the quantitative determination of the distribution of particle sizes in soils. The distribution of particle sizes larger than 75 μ m (retained on the No. 200 sieve) is determined by sieving, while the distribution of particle sizes smaller than 75 μ m is determined by a sedimentation process, using a hydrometer to secure the necessary data (Notes 1 and 2).

Note 1—Separation may be made on the No. 4 (4.75-mm). No. 40 (425- μ m), or No. 200 (75- μ m) sieve instead of the No. 10. For whatever sieve used, the size shall be indicated in the report.

Note 2—Two types of dispersion devices are provided: (1) a high-speed mechanical stirrer, and (2) air dispersion. Extensive investigations indicate that air-dispersion devices produce a more positive dispersion of plastic soils below the 20- μ m size and appreciably less degradation on all sizes when used with sandy soils. Because of the definite advantages favoring air dispersion, its use is recommended. The results from the two types of devices differences in particle size distribution, especially for sizes finer than 20 μ m.

2. Referenced Documents

2.1 ASTM Standards:

- 7 D421 Practice for Dry Preparation of Soil Samples for Particle-Size Analysis and Determination of Soil Constants²
 - E 11 Specification for Wire-Cloth Sieves for Testing Purposes³
 - E 100 Specification for ASTM Hydrometers⁴

3. Apparatus

3.1 Balances — A balance sensitive to 0.01 g for weighing the material passing a No. 10 (2.00-mm) sieve, and a balance sensitive to 0.1 % of the mass of the sample to be weighed for weighing the material retained on a No. 10 sieve.

3.2 Stirring Apparatus—Either apparatus A or B may be used.

3.2.1 Apparatus A shall consist of a mechanically operated stirring device in which a suitably mounted electric motor turns a vertical shaft at a speed of not less than 10 000 rpm without load. The shaft shall be equipped with a

* Annual Book of ASTM Standards, Vol 14.01.

replaceable stirring paddle made of metal. plastic. or hard rubber. as shown in Fig. 1. The shaft shall be of such length that the stirring paddle will operate not less than $\frac{1}{2}$ in. (19.0 mm) nor more than $\frac{1}{2}$ in. (38.1 mm) above the bottom of the dispersion cup. A special dispersion cup conforming to either of the designs shown in Fig. 2 shall be provided to hold the sample while it is being dispersed.

3.2.2 Apparatus B shall consist of an air-jet dispersion cup⁵ (Note 3) conforming to the general details shown in Fig. 3 (Notes 4 and 5).

Note 3—The amount of air required by an air-jet dispersion cup is of the order of 2 R^3 /min: some small air compressors are not capable of supplying sufficient air to operate a cup.

Note 4—Another air-type dispersion device, known as a dispersion tube, developed by Chu and Davidson at Iowa State College, has been shown to give results equivalent to those secured by the air-jet dispersion cups. When it is used, soaking of the sample can be done in the sedimentation cylinder, thus eliminating the need for transferring theslurry. When the air-dispersion tube is used, it shall be so indicated in the report.

NOTE 5—Water may condense in air lines when not in use. This water must be removed, either by using a water trap on the air line, or by blowing the water out of the line before using any of the air for dispersion purposes.

3.3 Hydrometer—An ASTM hydrometer. graduated to read in either specific gravity of the suspension or grams per litre of suspension, and conforming to the requirements for hydrometers 151H or 152H in Specifications E 100. Dimensions of both hydrometers are the same, the scale being the only item of difference.

3.4 Sedimentation Cylinder—A glass cylinder essentially 18 in. (457 mm) in height and 2½ in. (63.5 mm) in diameter, and marked for a volume of 1000 mL. The inside diameter shall be such that the 1000-mL mark is 36 ± 2 cm from the bottom on the inside.

3.5 Thermometer—A thermometer accurate to I'F (0.5°C).

3.6 Sieves—A series of sieves, of square-mesh woven-wire cloth, conforming to the requirements of Specification E 11. A full set of sieves includes the following (Note 6):

⁴ This method is under the jurisdiction of ASTM Committee D-18 on Soil and Rock and is the direct responsibility of Subcommittee D18.03 on Texture. Plasticity, and Density Characteristics of Soils.


Current edition approved Nov. 21, 1963. Originally published 1935. Replaces D 422 - 62.

² Annual Book of ASTM Standards, Vol 04.08.

³ Annual Book of ASTM Standards. Vol 14.02.

⁹ Detailed working drawings for this cup are available at a nominal cost frothe American Society for Testing and Materials, 1916 Race St., Philadelphia, P. 19103. Order Adjunct No. 12-404220-00.



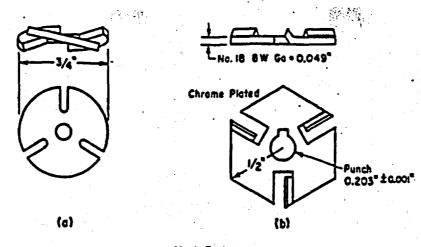

	Metric Equivalenta								
.	0 001	0.049	0.203	1/2	4.	_			
mm	0.03	1.24	5.16	12.7	19.0				

FIG. 1 Detail of Stirring Paddles

. 3-in. (75-mm) 2-in. (50-mm) 1*5-in. (37 5-mm) 1-in. (25.0-mm) 2*V⊷in. (19.0-mm) 3*™⊷in. (9.5-mm) 3*No. 4 (4.75-mm)

Note 6-A set of sieves giving uniform spacing of points for the graph. as required in Section 17, may be used if desired. This set consists of the following sieves:

(3-in. (75-mm) (34-in. (37.5-mm) 84-in. (19.0-mm) 7-1-in. (9.5-mm) (No. 4 (4.75-mm) (No. 8 (2.36-mm) No. 16 (1.18-mm) No. 30 (600-µm) No. 50 (300-µm) No. 100 (150-µm) No. 200 (75-µm)

No. 10 (2.00-mm)

No. 20 (850-um)

No. 40 (425-µm)

No. 60 (250-um)

No. 200 (75-µm)

No. 140 (106-µm)

3.7 Water Bath or Constant-Temperature Room—A water bath or constant-temperature room for maintaining he soil suspension at a constant temperature during the hydrometer analysis. A satisfactory water tank is an insulated tank that maintains the temperature of the suspension at a convenient constant temperature at or near 68°F (20°C). Such a device is illustrated in Fig. 4. In cases where the work is performed in a room at an automatically controlled constant temperature, the water bath is not necessary.

3.8 Beaker—A beaker of 250-mL capacity.

3.9 Timing Device—A watch or clock with a second hand.

J. Dispersing Agent

4.1 A solution of sodium hexametaphosphate (sometimes called sodium metaphosphate) shall be used in distilled or temineralized water, at the rate of 40 g of sodium texametaphosphate/litre of solution (Note 7).

NOTE 7—Solutions of this salt, if acidic, slowly revert or hydrolyze ack to the orthophosphate form with a resultant decrease in dispersive "n. Solutions should be prepared frequently (at least once a month) justed to pH of 8 or 9 by means of sodium carbonate. Bottles ng solutions should have the date of preparation marked on "etc..

4.2 All water used shall be either distilled or demineralized water. The water for a hydrometer test shall

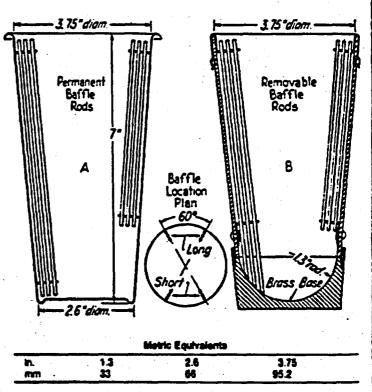


FIG. 2 Dispersion Cups of Apperatus

be brought to the temperature that is expected to prevail during the hydrometer test. For example, if the sedimentation cylinder is to be placed in the water bath, the distilled or demineralized water to be used shall be brought to the temperature of the controlled water bath; or, if the sedimentation cylinder is used in a room with controlled temperature, the water for the test shall be at the temperature of the room. The basic temperature for the hydrometer test is 68[°]F (20°C). Small variations of temperature do not introduce differences that are of practical significance and do not prevent the use of corrections derived as prescribed.

AR 301421

5. -

out

dur

con

siev

the

Dur

suff

folic

sha!

the

S.

S. be a;

g for

weig

Scog

Wash

Dere

Cala

Pulve

Dessin

and c

6. P

mm)

6.;

Nc

S.

5

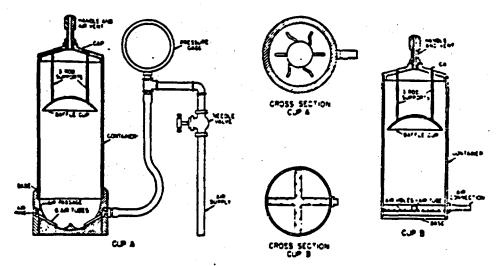


FIG. 3 Air-Jet Dispersion Cupe of Apparatus B

5. Test Sample

۴.

لتهر

nta

100

5.1 Prepare the test sample for mechanical analysis as outlined in Practice D 421. During the preparation procedure the sample is divided into two portions. One portion contains only particles retained on the No. 10 (2.00-mm) sieve while the other portion contains only particles passing the No. 10 sieve. The mass of air-dried soil selected for purpose of tests, as prescribed in Practice D 421, shall be sufficient to yield quantities for mechanical analysis as follows:

5.1.1 The size of the portion retained on the No. 10 sieve shall depend on the maximum size of particle, according to the following schedule:

Nominal Diameter of Largest Particles. ia. (mm)	Approximate Minimum Mass of Portuon, g
¥a (9.5)	500
¥• (19.0)	1000
1 (25.4)	2000
11/2 (38.1)	3000
2 (50.8)	4000
3 (76.2)	5000

5.1.2 The size of the portion passing the No. 10 sieve shall be approximately 115 g for sandy soils and approximately $\delta 5$ g for silt and clay soils.

5.2 Provision is made in Section 5 of Practice D 421 for weighing of the air-dry soil selected for purpose of tests, the separation of the soil on the No. 10 sieve by dry-sieving and washing, and the weighing of the washed and dried fraction retained on the No. 10 sieve. From these two masses the percentages retained and passing the No. 10 sieve can be calculated in accordance with 12.1.

Note 5—A check on the mass values and the thoroughness of pulverization of the clods may be secured by weighing the portion passing the No. 10 sieve and adding this value to the mass of the washed and oven-dried portion retained on the No. 10 sieve.

SIEVE ANALYSIS OF PORTION RETAINED ON NO. 10 (2.00-mm) SIEVE

÷.

6. Procedure

6.1 Separate the portion retained on the No. 10 (2.00mm) sieve into a series of fractions using the 3-in. (75-mm),

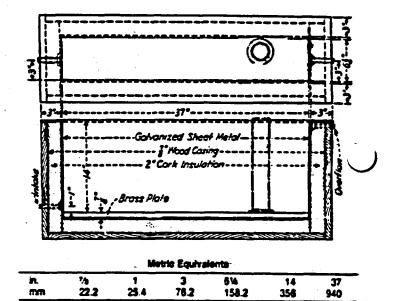


FIG. 4 Insulated Water Bath

2-in. (50-mm), 1½-in. (37.5-mm), 1-in. (25.0-mm), ¾-in. (19.0-mm), ¾-in. (9.5-mm), No. 4 (4.75-mm), and No. 10 sieves, or as many as may be needed depending on the sample, or upon the specifications for the material under test.

6.2 Conduct the sieving operation by means of a lateral and vertical motion of the sieve, accompanied by a jarring action in order to keep the sample moving continuously over the surface of the sieve. In no case turn or manipulate fragments in the sample through the sieve by hand. Continue sieving until not more than 1 mass % of the residue on a sieve passes that sieve during 1 min of sieving. When mechanical sieving is used, test the thoroughness of sieving by using the hand method of sieving as described above?

6.3 Determine the mass of each fraction on a balance conforming to the requirements of 3.1. At the end of weighing, the sum of the masses retained on all the sieves used should equal closely the original mass of the quantity sieved.

AR 301422

HYDROMETER AND SIEVE ANALYSIS OF PORTION PASSING THE NO. 10 (2.00-mm) SIEVE

Determination of Composite Correction for Hydrometer-Reading

7.1 Equations for percentages of soil remaining in suspension, as given in 14.3, are based on the use of distilled or demineralized water. A dispersing agent is used in the water, however, and the specific gravity of the resulting liquid is appreciably greater than that of distilled or demineralized water.

7.1.1 Both soil hydrometers are calibrated at 68°F (20°C). and variations in temperature from this standard temperature produce inaccuracies in the actual hydrometer readings. The amount of the inaccuracy increases as the variation from the standard temperature increases.

7.1.2 Hydrometers are graduated by the manufacturer to be read at the bottom of the meniscus formed by the liquid on the stem. Since it is not possible to secure readings of soil suspensions at the bottom of the meniscus, readings must be taken at the top and a correction applied.

7.1.3 The net amount of the corrections for the three items enumerated is designated as the composite correction, and may be determined experimentally.

7.2 For convenience, a graph or table of composite corrections for a series of 1° temperature differences for the range of expected test temperatures may be prepared and used as needed. Measurement of the composite corrections may be made at two temperatures spanning the range of expected test temperatures, and corrections for the intermee temperatures calculated assuming a straight-line rela-

Laship between the two observed values.

7.3 Prepare 1000 mL of liquid composed of distilled or demineralized water and dispersing agent in the same proportion as will prevail in the sedimentation (hydrometer) test. Place the liquid in a sedimentation cyclinder and the cylinder in the constant-temperature water bath, set for one of the two temperatures to be used. When the temperature of the liquid becomes constant, insert the hydrometer, and, after a short interval to permit the hydrometer to come to the temperature of the liquid, read the hydrometer at the top of the meniscus formed on the stem. For hydrometer 151H the composite correction is the difference between this reading and one; for hydrometer 152H it is the difference between the reading and zero. Bring the liquid and the hydrometer to the other temperature to be used, and secure the composite correction as before.

8. Hygroscopic Moisture

8.1 When the sample is weighed for the hydrometer test, weigh out an auxiliary portion of from 10 to 15 g in a small metal or glass container, dry the sample to a constant mass in an oven at 230 \pm 9°F (110 \pm 5°C), and weigh again. Record the masses.

9. Dispersion of Soil Sample

I When the soil is mostly of the clay and silt sizes, weigh o. a sample of air-dry soil of approximately 50 g. When the soil is mostly sand the sample should be approximately 100

9.2 Place the sample in the 250-mL beaker and cover with 125 mL of sodium hexametaphosphate solution (40 g/L). Stir until the soil is thoroughly wetted. Allow to soak for at least 16 h.

9.3 At the end of the soaking period, disperse the sample further, using either stirring apparatus A or B. If stirring apparatus A is used, transfer the soil - water slurry from the beaker into the special dispersion cup shown in Fig. 2. washing any residue from the beaker into the cup with distilled or demineralized water (Note 9). Add distilled or demineralized water, if necessary, so that the cup is more than half full. Stir for a period of 1 min.

NOTE 9-A large size syringe is a convenient device for handling the water in the washing operation. Other devices include the wash-water bottle and a hose with nozzle connected to a pressurized distilled water tank.

9.4 If stirring apparatus B (Fig. 3) is used, remove the cover cap and connect the cup to a compressed air supply by means of a rubber hose. A air gage must be on the line between the cup and the control valve. Open the control valve so that the gage indicates 1 psi (7 kPa) pressure (Note 10). Transfer the soil - water slurry from the beaker to the air-jet dispersion cup by washing with distilled or demineralized water. Add distilled or demineralized water, if necessary, so that the total volume in the cup is 250 mL, but no more.

Note 10-The initial air pressure of 1 psi is required to prevent the soil - water mixture from entering the air-jet chamber when the mixture is transferred to the dispersion cup.

9.5 Place the cover cap on the cup and open the air control valve until the gage pressure is 20 psi (140 kPa). Disperse the soil according to the following schedule:

Plasticity Index	
	616
Under S	5
6 to 20	10
Over 20	15

Soils containing large percentages of mica need be dispersed for only 1 min. After the dispersion period, reduce the gage pressure to 1 psi preparatory to transfer of soil - water slurry to the sedimentation cylinder.

10. Hydrometer Test

10.1 Immediately after dispersion, transfer the soil - water slurry to the glass sedimentation cylinder, and add distilled or demineralized water until the total volume is 1000 mL.

10.2 Using the palm of the hand over the open end of the cylinder (or a rubber stopper in the open end), turn the cylinder upside down and back for a period of 1 min 10 complete the agitation of the slurry (Note 11). At the end of 1 min set the cylinder in a convenient location and take hydrometer readings at the following intervals of time (measured from the beginning of sedimentation), or as many as may be needed, depending on the sample or the specification for the material under test: 2, 5, 15, 30, 60, 250, and 1440 min. If the controlled water bath is used, the sediment tation cylinder should be placed in the bath between the 2and 5-min readings.

Note 11-The number of turns during this minute should be approximately 60, counting the turn upside down and back as two turns

DOL: for und 12. : 12 divid origiı by 10 the m 12. (4.75. No. 1 and n soil p: SOU D: the 3/2 Remai: manni 12.3 Siève. mass c 13. H 13.1 ratio b

lir-dry

except

14. Per

14.1

11

10

¢î,

řτ.

th

12

in

1Å

for

100

SUS

\$10

11.

the

wat

the

230

1

hydron: hygrosc AR 301423

Any soil remaining in the bottom of the cylinder during the first few turns should be loosened by vigorous shaking of the cylinder while it is in the inverted position.

10.3 When it is desired to take a hydrometer reading, carefully insert the hydrometer about 20 to 25 s before the reading is due to approximately the depth it will have when the reading is taken. As soon as the reading is taken, carefully remove the hydrometer and place it with a spinning motion in a graduate of clean distilled or demineralized water.

Note 12—It is important to remove the hydrometer immediately after each reading. Readings shall be taken at the top of the meniscus formed by the suspension around the stem, since it is not possible to secure readings at the bottom of the meniscus.

10.4 After each reading, take the temperature of the suspension by inserting the thermometer into the suspension-

11. Sieve Analysis

. ∕le

1.72

i the

: 2.

nore

i he

water

water

"i une

y by

I ae

rj ol

Note

). the

1 or

ð., if

, but

n be

xture

____d

gage

i ter

timed

) he

in to

nd of

nany

io de turns

° *

nd

ke L.ne

i he

nL.

TY

⁺th ∟or

> 11.1 After taking the final hydrometer reading, transfer the suspension to a No. 200 (75- μ m) sieve and wash with tap water until the wash water is clear. Transfer the material on the No. 200 sieve to a suitable container, dry in an oven at 230 \pm 9°F (110 \pm 5°C) and make a sieve analysis of the portion retained, using as many sieves as desired, or required for the material. or upon the specification of the material under test.

CALCULATIONS AND REPORT

12. Sieve Analysis Values for the Portion Coarser than the No. 10 (2.00-mm) Sieve

12.1 Calculate the percentage passing the No. 10 sieve by dividing the mass passing the No. 10 sieve by the mass of soil originally split on the No. 10 sieve, and multiplying the result by 100. To obtain the mass passing the No. 10 sieve, subtract the mass retained on the No. 10 sieve from the original mass.

12.2 To secure the total mass of soil passing the No. 4 (4.75-mm) sieve, add to the mass of the material passing the No. 10 sieve the mass of the fraction passing the No. 4 sieve and retained on the No. 10 sieve. To secure the total mass of soil passing the %-in. (9.5-mm) sieve, add to the total mass of soil passing the No. 4 sieve, the mass of the fraction passing the %-in. sieve and retained on the No. 4 sieve. For the remaining sieves, continue the calculations in the same manner.

12.3 To determine the total percentage passing for each sieve, divide the total mass passing (see 12.2) by the total mass of sample and multiply the result by 100.

13. Hygroscopic Moisture Correction Factor

13.1 The hydroscopic moisture correction factor is the ratio between the mass of the oven-dried sample and the air-dry mass before drying. It is a number less than one, except when there is no hygroscopic moisture.

14. Percentages of Soil in Suspension

14.1 Calculate the oven-dry mass of soil used in the hydrometer analysis by multiplying the air-dry mass by the hygroscopic moisture correction factor.

TABLE	Values of Correction Factor, a, for Different Specific
	Gravities of Soil Particles ⁴

Specific Gravity	Conscion Factor*
2.95	0.94
2.90	0.95
2.85	0.96
2.80	0.97
2.75	0.98
2.70	0 99
2.65	1.00
2.50	1.01
2.55	1 02
2.50	1.03
2 45	1.05

⁴ For use in equilition for percentage of soil remaining in suspension when using Hydrometer 152H.

14.2 Calculate the mass of a total sample represented by the mass of soil used in the hydrometer test, by dividing the oven-dry mass used by the percentage passing the No. 10 (2.00-mm) sieve, and multiplying the result by 100. This value is the weight W in the equation for percentage remaining in suspension.

14.3 The percentage of soil remaining in suspension at the level at which the hydrometer is measuring the density of the suspension may be calculated as follows (Note 13): For hydrometer 151H:

$P = [(100\ 000/\mathcal{W}) \times G/(G - G_1)](R - G_1)$

Note 13—The bracketed portion of the equation for hydrometer 151H is constant for a series of readings and may be calculated first and then multiplied by the portion in the parentheses.

where:

For hydrometer 152H:

- $P = (Ra/W) \times 100$
- a = correction faction to be applied to the reading of hydrometer 152H. (Values shown on the scale are computed using a specific gravity of 2.65. Correction factors are given in Table 1),
- P = percentage of soil remaining in suspension at the level at which the hydrometer measures the density of the suspension,
- R = hydrometer reading with composite correction applied (Section 7),
- W = oven-dry mass of soil in a total test sample represented by mass of soil dispersed (see 14.2), g,
- G = specific gravity of the soil particles, and
- G_i = specific gravity of the liquid in which soil particles are suspended. Use numerical value of one in both instances in the equation. In the first instance any possible variation produces no significant effect, and in the second instance, the composite correction for R is based on a value of one for G_i .

15. Diameter of Soil Particles

15.1 The diameter of a particle corresponding to the percentage indicated by a given hydrometer reading shall be calculated according to Stokes' law (Note 14), on the basis that a particle of this diameter and the supersion at the beginning of sedimentation and had settler to the level at which the hydrometer is measuring the density of the suspension. According to Stokes' law:

 $D = \sqrt{130\pi/980(G - G.11 \times L/T)}$

TABLE 2

where:

- D = diameter of particle, mm,
 - coefficient of viscosity of the suspending medium (in this case water) in poises (varies with changes in temperature of the suspending medium),
- L = distance from the surface of the suspension to the level at which the density of the suspension is being measured, cm. (For a given hydrometer and sedimentation cylinder, values vary according to the hydrometer readings. This distance is known as effective depth (Table 2)).
- T = interval of time from beginning of sedimentation to the taking of the reading, min,
- G = specific gravity of soil particles, and
- G₁ = specific gravity (relative density) of suspending medium (value may be used as 1.000 for all practical purposes).

NOTE 14—Since Stokes' law considers the terminal velocity of a single sphere falling in an infinity of liquid, the sizes calculated represent the diameter of spheres that would fall at the same rate as the soil particles.

15.2 For convenience in calculations the above equation may be written as follows:

 $D = K \sqrt{L/T}$

where:

K = constant depending on the temperature of the suspension and the specific gravity of the soil particles. Values of K for a range of temperatures and specific gravities are given in Table 3. The value of K does not change for a series of readings constituting a test, while values of L and T do vary.

15.3 Values of D may be computed with sufficient accuracy, using an ordinary 10-in. slide rule.

NOTE 15—The value of L is divided by T using the A- and B-scales, the square root being indicated on the D-scale. Without ascertaining the value of the square root it may be multiplied by K, using either the C- or CI-scale.

Sieve Analysis Values for Portion Finer than No. 10 (2.00-mm) Sieve

16.1 Calculation of percentages passing the various sieves used in sieving the portion of the sample from the hydrometer test involves several steps. The first step is to calculate the mass of the fraction that would have been retained on the No. 10 sieve had it not been removed. This mass is equal to the total percentage retained on the No. 10 sieve (100 minus total percentage passing) times the mass of the total sample represented by the mass of soil used (as calculated in 14.2), and the result divided by 100.

16.2 Calculate next the total mass passing the No. 200 sieve. Add together the fractional masses retained on all the sieves, including the No. 10 sieve, and subtract this sum from the mass of the total sample (as calculated in 14.2).

16.3 Calculate next the total masses passing each of the

r sieves, in a manner similar to that given in 12.2. .6.4 Calculate last the total percentages passing by dividing the total mass passing (as calculated in 16.3) by the total mass of sample (as calculated in 14.2), and multiply the result by 100.

Hydromet	# 151H	Hydrometer 152H				
Actual	Effective		Effective	Actual	Effective	
Hydrometar -	Depth.	Hydrometer	Depth.	Hydrometer	Depth.	
Reading	L. cm	Reading	L, ch	Reading	L. em	
1.000	16.3	a	16.3	31	11.2	
1.001	16.0	1	16.1	32	11.1	
1.002	15.8	2	16.0	33	109	
1.003	15.5	3	158	34	10.7	
1.004	15.2	4	15.6	35	10.6	
1 005	15.0	\$	15 5	•		
1.005	14.7	Ē	15.3	36	10.4	
1.007	14.4	7	15.2	37	10.2	
1.008	14.2		15.0	38	10.1	
1.009	13.9	.	14 8	39	9.9	
1.010	13.7	10	14.7	40	9.7	
1.011	13.4	11		. 41		
1.012	13.4	11	14 S 14.3	41	9.6	
1.013	12.9	13	14.3	42	9.4 9.2	
1.014	12.6	. 14	14 0	4		
1.015	12.3	15	13.8	45	9.1 8.9	
				-		
1.016	12.1	16	13.7	48	8.8	
1.017	11.8	17	13.5	47	8.6	
1.018	115	18	13.3	48	8.4	
1.019	11.3	19	13.2	49	8.3	
1.020	11.0	20	13.0	50	8.1	
1.021	10.7	. 21	12.9	51	7.9	
1.022	10.5	22	12.7	52	7.8	
1.023	10.2	23	12.5	53	7.6	
1.024	10.0	24	12.4	54	7.4	
1.025	9.7	25	12.2	55	7.3	
1.025	9.4	26	12.0	56	7.1	
1.027	1.2	27	11.9	57	7.0	
1.028	8.9	28	11.7	58	6.6	
1.029	8.6	29	11.5	59	6.6	
1.030	6.4	30	11.4	60	6.5	
1.031	6.1					
1.032	7.8					
1.033	7.6			-		
1.034	7.3					
1.035	7.0					
1.036	6.8					
1.037	6.5					
1.038	6.2				1	
		opth are calculate		40 (8000)		
	-	L = L, + 1/2 [Lg				

ſ

£

F

2

2

0 fi

1.

٩.

()

2n

Nc Arr

dis

512

RE

6120

Spc.

sha.

shal

SDer

1;

I

.

Values of Effective Depth Based on Hydrometer and

Sedimentation Cylinder of Specified Sizes⁴

- Le · overall length of the hydrometer bub, on,
- $V_{\rm e}$ = volume of hydrometer bub, cm³, and $A_{\rm e}$ = cross-sectional area of sedimentation cylinder, cm²
- Values used in calculating the values in Table 2 are as follow
- For both hydrometers, 151H and 152H:
- L. = 14.0 cm
- Ve = 67.0 cm
- A = 27.8 cm² For hydrometer 151H:
- L. = 10.5 cm for a meding of 1.000
- = 2.1 on for a reading of 1.031 AR301425
- For hydrometer 152H: $L_1 = 10.5$ cm for a reading of 0 g/fore
- = 2.3 cm for a reading of 50 g/itre

17. Graph

17.1 When the hydrometer analysis is performed, a graph

[] D 422

Temperature.ar	Specific Gravity of Soll Particles										
•C	2.45	2.50	2.55	2.50	2.65	2.70	2.75	2.80	2.85		
15	. 0.01510	0.01505	0.01481	0.01457	0.01435	0.01414	0.01394	0.01374	0 01358		
17	0.01511	0.01486	0.01462	0.01439	0.01417	0.01396	0.01375	0.01358	0.01338		
18	0.01492	0.01457	0.01443	0.01421	0.01399	0.01378	0.01359	0.01339	0.01321		
19	0.01474	0.01449	0.01425	0.01403	0.01382	0.01361	0.01342	0.1323	0.01305		
20	0.01458	0.01431	0.01408	0.01386	0.01365	0.01344	0.01325	0.01307	0.01289		
21	0.01438	0.01414	0.01391	0.01369	0.01348	0.01328	0.01309	0.01291	0.01273		
22	0.01421	0.01397	0.01374	0.01353	0.01332	0.01312	0.01294	0.01276	0.01258		
23	0.01404	0.01381	0.01358	0.01337	0.01317	0.01297	0.01279	0.01261	0.01243		
24	0.01388	0.01365	0.01342	0.01321	0.01301	0.01282	0.01264	0.01246	0.01229		
24 25	0.01372	0.01349	0.01327	0.01305	0.01255	0.01257	0.01249	0.01232	0.01215		
25	0.01357	0.01334	0.01312	0.01291	0.01272	0.01253	0.01235	0.01218	0.01201		
25 27	0 01342	0.01319	. 0.01297	0.01277	0.01258	0.01239	0.01221	0.01204	0.01168		
28	0.01327	0.01304	0.01253	0.01254	0.01244	0.01255	0.01208	0.01191	0.01175		
29	0.01312	0.01290	0.01269	0.01249	0.01230	0.01212	0.01195	0.01178	0.01162		
30	0.01298	0.01278	0.01258	0.01236	0.01217	0.01199	0.01182	0.01165	0.01149		

TABLE 3 Values of K for Use in Equation for Computing Diameter of Particle in Hydrometer Analysia

of the test results shall be made, plotting the diameters of the particles on a logarithmic scale as the abscissa and the percentages smaller than the corresponding diameters to an arithmetic scale as the ordinate. When the hydrometer analysis is not made on a portion of the soil, the preparation of the graph is optional, since values may be secured directly from tabulated data.

18. Report

17

1

4

to the

18.1 The report shall include the following:

18.1.1 Maximum size of particles.

18.1.2 Percentage passing (or retained on) each sieve, which may be tabulated or presented by plotting on a graph (Note 16).

18.1.3 Description of sand and gravel particles:

18.1.3.1 Shape—rounded or angular,

18.1.3.2 Hardness-hard and durable, soft, or weathered and friable,

18.1.4 Specific gravity, if unusually high or low,

18.1.5 Any difficulty in dispersing the fraction passing the No. 10 (2.00-mm) sieve, indicating any change in type and amount of dispersing agent, and

18.1.6 The dispersion device used and the length of the dispersion period.

Note 16—This tabulation of graph represents the gradation of the sample tested. If particles larger than those contained in the sample were removed before testing, the report shall so state giving the amount and maximum size.

18.2 For materials tested for compliance with definite specifications, the fractions called for in such specifications shall be reported. The fractions smaller than the No. 10 sieve shall be read from the graph.

18.3 For materials for which compliance with definite specifications is not indicated and when the soil is composed

almost entirely of particles passing the No. 4 (4.75-mm) sieve, the results read from the graph may be reported as follows:

- (1) Gravel, passing 3-in. and retained on No. 4 sieve
- (2) Sand, passing No. 4 sieve and retained on No. 200 sieve
 (a) Coarse and, passing No. 4 sieve and retained on No. 10 sieve
- (3) Silt size, 0.074 to 0.005 mm
- (4) Clay size, smaller than 0.005 mm
- Colloids, smaller than 0.001 mm

18.4 For materials for which compliance with definite specifications is not indicated and when the soil contains material retained on the No. 4 sieve sufficient to require a sieve analysis on that portion, the results may be reported as follows (Note 17):

STEVE ANALYSIS-

Sieve Size-		Percentage Passing
3-in. 2-in.		• • • • • • • • • •
1%-in. 1-in. %-in.		················
%-in. No. 4 (4.75-mm)		•••••
No. 10 (2.00-mm) No. 40 (425-µm)		·• ·•• · · · · · ·
No. 200 (75-µm)	HYDROMETER ANALYSIS	••••
0.074 mm 0.005 mm 0.001 mm		••••••

NOTE 17-No. 8 (2.36-mm) and No. 50 (300-µm) sieves may be substituted for No. 10 and No. 40 sieves.

The American Society for Testing and Meterials takes no position respecting the validity of any patent rights asserted in connection with any term mentioned in this standard. Users of this standard are expressly advised that determination of the validity R any BCT 426 petert rights, and the risk of intringement of such rights, are antirely that own responsibility.

This standard is subject to revision at any time by the responsible technical committee and must be reviewed every five years and if not revised, either responsed or withdrawn. Your comments are invited either for revision of this standards and should be addressed to ASTM Headquarters. Your comments will receive careful consideration at a meeting of the responsible technical committee, which you may attend. If you feel that your comments have not received a fair hearing you should make your views known to the ASTM Committee on Standards, 1918 Rece St., Philadelphia, PA 19103.

- CARBON, TOTAL ORGANIC AND INCRGANIC

يذور والمتعاط

Carbon may exist in sediment and water samples as either inorganic or organic compounds. Inorganic carbon is present as carbonates, bicarbonates, and possibly free carbon dioxide. Specific types of compounds that are considered to be included in the organic carbon fraction are nonvolatile organic compounds (sugars), volatile organic compounds (mercaptans), partially volatile compounds (oils), and particulate carbonaceous materials (cellulose).^{1,2[#]}

The basis of the method is the catalytic or chemical oxidation of carbon in carbon-containing compounds to carbon dioxide followed by the quantification of the carbon dioxide produced. Alternately, the carbon may be reduced to methane and appropriately quantified. It follows, then, that the distinction between inorganic carbon and organic carbon is the method of sample pretreatment. There are presently two procedures for defining this separation. One method is based on sample treatment with a strong acid. Analysis of an untreated sample is a measure of total carbon while analysis of the acid-treated fraction is a measure of organic carbon. Inorganic carbon is based on differential thermal combustion with organic compounds being converted to carbon dioxide at $950^{\circ}C$ to $1300^{\circ}C$.^{4,5}

Sample Handling and Storage

Flowcharts for the handling of samples intended for organic carbon and inorganic carbon analysis are presented in Figure 3-6 and Figure 3-7. Water and sediment samples to be analyzed for inorganic carbon may be stored in glass or plastic containers. There is no effective preservative because of the carbon dioxide reserve in the atmosphere. The only precaution that can be taken for inorganic

* References for this procedure can be found on page 3-76. AR301427

Handling and storage of sumples for total inorganic curbon analysis 3 CORE SECTION Total Sadiment Conc. 95 ļ **6.** 4°C · CONE SAUNE DICESI 1015) 1015) 510 3 REDGE SAMPLE Nabile Conc. STORE LET ļ 21 12 **SIA** 18 19 ••• lised in Elutriate ļ 5 **.**. 4°C 4°C 4°C 4°C 4°C 4°C ELUTAIATE ANALYZE (51A) Filler 7 8 9 24 12 9 • ND TREATMENT 5 Total Vatar Conc. 78 99 1 24 PC 3 **4**°9 WIER SAMLE FILTER ALIGIA STORE (27) (27) Figure 3-6. SAMPLE VOLUME ON LEIGHT SUPLE DESIGNATION DICESTION SOLUTION SAMPLE TREATMENT PALEGRAATIVE ACIDIFY STORAGE TINE DICEST STORE (11) (11) CONTAINER PUMPOSE •

3

ļ

AR301428

CONE SAMLE (000 STONE (000 (000) (\$10 \$2 \$3	Total Total Total Sediment Sediment Conc. Conc. Conc.	£." 9 £."	None Air dry Freeze	None None		oo el 30-100 el Vertable 0.5-5.0 g 0.5-5.0 g 0.5-5.0 g 0.5-5.0 g Handling and storage of samples for total organic carbon analysis
	SIA	Mobi la Conc.	4. U	Rone	7 4	ŝ	uf samples fo
	5	tsed in Elutriate	u	ł	4.C	2	vertable OTAGE O
	5	Soluble Nater Cenc.	.	filter	IKI or N ₂ 50, to pH2 4°C	24 hr	æ-tœ =t ing and st
	5	Total Nator Conc.	پ		NC1 or N ₁ 50, to pH2 405	a z	Ē
	SWPLE DESIGNATION	juanose	CONTAINER	SAMPLE TREATMENT	PRESERVATIVE	AR30	Figure 3-7.

44日十

carbon is to completely fill the sample container at the time of sampling (exclude all air bubbles), tightly seal the container, and complete the analysis immediately (Figure 3-6).

Water samples for organic carbon analysis should be stored in glass containers unless substitute containers have been shown not to affect total organic carbon (TOC) analyses. Samples should be processed as soon as possible (within 24 hr if possible) to minimize change due to chemical or biological oxidation. Atmospheric uptake of carbon dioxide is less critical since it would be evolved when the sample is acidified prior to analysis. Sediment samples for organic carbon analysis may be stored in either plastic or glass containers (Figure 3-7). Air drying of sediments (S2) may lead to low TOC results due to oxidation or volatilization. Therefore, moist storage (SID) or frozen storage (S3) would be the preferred method of storage. If samples are frozen, excessive temperatures should not be used to thaw the samples.

3-68

Procedure for Water Samples (W1, W2, S1A)

Method 1: Infrared Analysis 6,7

Apparatus

Sample homogenizer such as a Waring blender or ultrasonic blender Magnetic stirrer

Hypodermic syringe

Total carbon analyzer, either a single channel or a dual channel instrument (Dow-Beckman Carbonaceous Analyzer Model No. 915, Dohrmann Envirotech DC-50 carbon analyzer, Oceanography International Total Carbon Analyzer, Leco, or equivalent)

Reagents

Distilled water: the distilled water used in the preparation of standards and dilution of samples should be of the highest quality in order to have a small blank.

Organic carbon, stock solution, 1000 mg/f C: dissolve 2.125 g anhydrous potassium biphthalate, KHC: H:O:, in distilled water and dilute to 1 ^f in a volumetric flask.

Organic carbon, standard solutions: prepare standard solution by dilution of the stock solution as required.

Inorganic carbon, stock solution, 1000 mg/f: dissolve 3.500 g sodium bicarbonate, NaHCO3, and 4.418 g sodium carbonate, Na2CO3, in distilled water in a 1-2 volumetric flask and make up to the mark.

Inorganic carbon, standard solution: prepare standards from the stock solution as required.

Facking for total carbon tube: dissolve 20 g cobalt nitrate, $CO(NO_3)_2$ · GH_2O_3 in 50 ml distilled water. Add this solution to 15 g long-fiber asbestos in a porcelain evaporating dish. Mix and evaporate to dryness on a steam bath. Place the dish in a muffle furnace and bring to 950°C. After 1 to 2 hr at this temperature, remove the dish and allow to cool. Break up any large lumps and mix adequately but not excessively. With the combustion tube held in a vertical position, taper joint up, put about 1/2 in. of untreated asbestos in the tube first, then transfer in small amounts, approximately 1 g of catalyst into the tube with forceps or tweezers. As it is added, tap or push the material gently with a 1/4-in. glass rod. Do not force the packing. The weight of the rod itself is sufficient to compress the material. When completed, the length of the packing should be about 5 or 6 cm. Test the packed tube by measuring the flow rate [3] [43]

ست و و و

gas through it at room temperature, and then at 750°C. The rate should not drop more than 20 percent.

Packing for carbonate tube (dual channel instrument): place a small wad of quartz wool or asbestos near the exit end of the carbonate evolution tube. From the entrance end add 6 to 12 mesh quartz chips, allowing these to collect against the wad to a length of 10 cm. Pour an excess of 85 percent phosphoric acid, H₃PO₄, into the tube while holding it vertically and allow the excess to drain out.

Nitrogen gas, carbon dioxide free.

Procedure

Turn on the infrared analyzer, recorder, and tube furnaces, setting the total carbon furnace at 950°C and the carbonate furnace at 175°C. Allow sufficient warm-up time for stable, drift-free operation; about 2 hr is required. If used daily, the analyzer can be left on continuously. Adjust the oxygen flow rate to 80 to 100 ml/min through the total carbon tube. With other instruments, follow manufacturer's directions to warm up the instrument.

Immediately prior to carrying out calibrations or analyses, inject several portions of the appropriate standard into the tube to be used, until constant readings are obtained. The actual injection technique is as follows: rinse the syringe several times with the solution to be analyzed, fill, and adjust the volume to be pipeted. Wipe off the excess with soft paper tissue, taking care that no lint adheres to the needle. Remove the plug from the syringe holder, insert the sample syringe, and inject the sample into the combustion tube with a single, rapid movement of the thumb. Leave the syringe in the holder until the flow rate returns to normal, then replace it with the plug.

Successively introduce a convenient sized aliquot (20 to $50 \ \mu$) of each organic carbon standard and a blank into the total carbon tube and record peak heights. Between injections allow the recorder pen to return to its baseline. When a dual channel instrument is used, the standardization procedure must be repeated using carbonate standards to calibrate the low temperature channel.

Thoroughly mix the sample. Inject a convenient sized aliquot (20 to 50 µl) of the sample into the total carbon tube and AR301432

record the peak height. This result is a measure of the organic carton concentration and the inorganic carbon concentration of the sample.

Thoroughly mix the sample using a Waring blender or an ultrasonic homogenizer. Transfer 10 to 15 ml of sample to a 30-ml beaker and acidify with concentrated HCl to a pH of 2 or less. Purge the sample with carbon dioxide free nitrogen gas for 5 to 10 min. Plastic tubing should not be used during the purging process unless it has been previously shown that it will not add organic carbon to the sample.

Mix the acidified sample on a magnetic stirrer. While stirring, withdraw a subsample from the beaker using a hypodermic needle with a 150-µm opening. Inject the sample into the carbon analyzer to be used and record the peak height. This result is a measure of the organic carbon concentration of the sample.

Using either clear or filtered water samples, analytical precision will approach 1 to 2 percent or 1 to 2 mg/l carbon, whichever is greater. Analytical precision for unfiltered water samples will increase to 5 to 10 percent because of the difficulty associated with sampling particulate matter and the fact that the needle opening of the syringe limits the maximum size of the particles that can be included in the sample.

Calculations

<u>Dual-channel instrument.</u> Prepare calibration curves derived from the peak heights obtained with the standard total carbon and inorganic carbon solutions.

Determine the concentration of total carbon and inorganic carbon in the sample by comparing sample peak heights with the calibration curves.

Determine the concentration of total inorganic carbon in the sample by subtracting the organic carbon value from the total carbon value.

Single-channel instrument. Prepare a calibration curve derived from the peak heights obtained with the standard total carbon solutions. Determine the total carbon concentration in the sample by comparing the peak height of the first sample injection with $\beta = 30 + 433$

calibration curve. Determine the organic carbon concentration in the sample by comparing the peak height of the second sample injection with the calibration curve. Inorganic carbon concentrations are calculated by subtracting the organic carbon concentration from the total carbon concentration.

K

Procedures for Sediment Samples (S1D, S3)

Method 1: Sample Ignition

Apparatus

Induction furnace such as the Leco WR-12, Dohrmann DC-50, Coleman CH analyzer, or Perkin Elmer 240 elemental analyzer

Combustion boats

Microbalance

Desiccator

Reagents

10 percent hydrochloric acid: mix 100 ml concentrated HCl with 900 ml distilled water.

Copper oxide fines.

Benzoic acid.

Procedure

Dry at 70°C and grind the sediment sample.

Weigh a combustion boat and record the weight. Place 0.2 to 0.5 g homogenized sediment in the combustion boat and reweigh. Combustion boats should not be handled with the bare hand during this process.

If total carbon or inorganic carbon is to be determined, Cupric oxide fines may be added to the sample to assist in combustion. Combust the sample in an induction furnace. Record the result as total carbon.

If organic carbon is to be determined, treat a known weight of dried sediment with several drops of 10 percent HCL. Wait until the effervescing is completed and add more acid. Continue this process until the incremental addition of acid causes no further effervescence. Do not add too much acid at one time as this may cause loss of sample due to frothing.

Dry the sample at 70° C and place in a desiccator. Add Cupric oxide fines, combust the sample in an induction furnace, and record the result as organic carbon.

AR301435

•

Calculations.

The carbon content of the sample can be calculated as:

%C = weight of tube (after-before) x 27.29 sample weight

Derivation of factor:

27.29 = 12.011 (molecular weight carbon) 44.011 (molecular weight carbon dioxide) * 100%

When the total sample results are used, the result is percent carbon in the sample. When acid-treated samples are used, the result is percent organic carbon. Inorganic carbon is calculated as total carbon minus organic carbon.

Method 2: Differential Combustion***

Apparatus

.....

.....

Sargent programmed microcombustion apparatus or equivalent Microbalance

Procedure

Air dry the sediment sample. Using a mortar and pestle, grind the sample to pass a 100-mesh screen.

Combust a known weight of sediment at a programmed heating rate of 300° to 950° C in 10 min and then maintain 950° C for 20 min. Trap the CO₂ in ascarite and record the weight as total carbon. A sample size should be selected that will produce 25 to 50 mg CO₂.

Weigh a second portion of the dried sediment. Combust this sample at a programmed rate of 300° to 650° C in 10 min and maintain 650° C for 20 min. Trap the CO₂ in ascarite and record the weight as organic carbon.

Calculations

The total carbon concentration, C_t , of the sample (in mg/g) is calculated as follows:

$$C_t = \frac{(x_t) (\frac{12}{55})}{(g)}$$

where

 x_t = weight of CO₂ evolved at 950°C, mg

3-74

g = weight of sample combusted, g

The organic carbon, C_0 , concentration of the sample (in $\pm \frac{\pi}{5}$) is calculated as follows:

.

 $C_{0} = \frac{(X_{0})(\frac{12}{64})}{(g)}$

where

.

 x_{o} = weight of CO₂ evolved at 650°C, mg

g = weight of sample combusted, g

Inorganic carbon, C_{I} , (in mg/g) is calculated as: $C_{I} = C_{t} - C_{o}$

Method 3: Wet Combustion 4.8

A third method has been used for carbon in sediments. This is based on the oxidation of the sample with dichromate and back titration of the sample with ferrous ammonium sulfate. References are provided for the procedure but details are not given. The procedure is similar to the chemical oxygen demand test which is not specific for carbon. The wet combustion method is a redox procedure and any reduced chemicals in the sediment samples (ferrous iron, manganous manganese, sulfide) will react with the dichromate. Therefore, this procedure is not recommended unless other instrumentation is not available.

3-75

References

- 1. U. S. Environmental Protection Agency. "Manual of Methods for Chemical Analysis of Water and Wastes." Methods Development and Quality Assurance Research Laboratory, National Environmental Research Center; Cincinnati, Ohio. 298 p. (1974).
- 2. U. S. Environmental Protection Agency. "Methods for Chemical Analysis of Water and Wastes." Environmental Monitoring and Support Laboratory, Office of Research and Development, EPA; Cincinnati, Ohio (1979).
- 3. Giovannini, G., Poggio, G., and Sequi, P. "Use of an Automatic CHN Analyzer to Determine Organic and Inorganic Carbon in Soils." Unpublished Report, Laboratory of Soil Chemistry, via Corridoni, Pisa, Italy. 9 p. (1975).
- Konrad, J. G., Chesters, G., and Keeney, D. R. "Determination of Organic- and Carbonate-Carbon in Freshwater Lake Sediments by a Microcombustion Procedure." J. Thermal Analysis 2:199-208 (1970).
- 5. Kemp, A. L. W. "Organic Matter in the Sediments of Lakes Ontario and Erie." Proc. 12th Conference Great Lakes Research 12:237-249 (1969).
- 6. Environment Canada. "Analytical Methods Manual," Inland Waters Directorate, Water Quality Branch; Ottawa, Canada (1974).
- American Public Health Association. <u>Standard Methods for the</u> <u>Examination of Water and Wastewater</u>. APHA; New York, New York. 1193 p. (1976).
- Gaudette, H. E., Flight, W. R., Toner, L., and Folger, D. W. "An Inexpensive Titration Method for the Determination of Organic Carbon in Recent Sediments." J. Sed. Petrology 44:249-253 (1974).